


CURAPROX

HAPPY PATIENT, HAPPY DENTIST

The greatest medicine of
all is teaching people how
not to need it.

Hippocrates

SWISS PREMIUM 
ORAL CARE

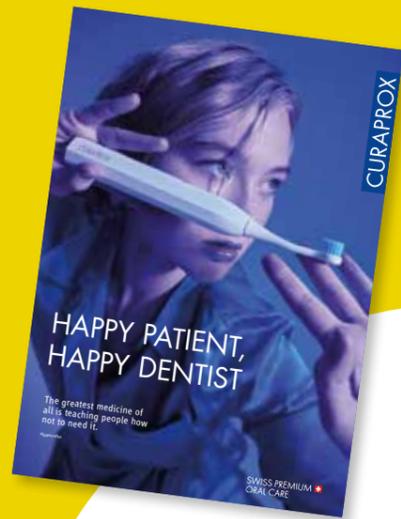

THIS CATALOGUE LEAVES (NOT) A MARK

We care about oral health. And we do so about our environment. That's why paper, inks and additives of this catalogue are fully recyclable and absolutely non-toxic. In accordance with the 'cradle to cradle®' approach, they can be completely returned to their corresponding material cycles at the end of their 'life'. The result: 0% waste - and this over several product cycles. That's the big difference between 'cradle to cradle®' and conventional recycling. A clean job, isn't it?



Höchster Standard für Ökoeffektivität.
Cradle to Cradle Certified™-Druckprodukte
hergestellt durch die Vögel AG.
Bindung ausgenommen.

gedruckt in der
schweiz



WORKING TOGETHER FOR THE BEST ORAL HEALTH – THANK YOU FOR OUR VALUABLE PARTNERSHIP

Good health is priceless. Too often, we do take a healthy, perfectly functioning body for granted. This has been brought into sharp focus during this Coronavirus pandemic.

As a dental professional, you can make a valuable contribution to ensuring that a healthy body remains healthy and is not burdened by issues that can be minimized or completely removed by appropriate means and the right methods.

The Curaden family and you, as a dental professional, share the same focus and passion. People, their oral health and their overall health are the centre of our world. Together we make an exceptional team: Curaden creates biofunctional, effective products that quickly bring visible and measurable success – and that are so much fun to use. You, meanwhile, deal directly with the patients, taking care of their oral health and making sure they use the products in the right way. This is the winning team in the fight against caries, bacteria and chronic diseases.

We are constantly delighted by your contribution and we will strive to do everything we can to make our partnership even more successful. We are obsessed with developing the very best products – thus offering you the ultimate support in your role as a dental professional, teacher and coach.

It is thanks to our tireless product development that we've discovered 'Perio plus'. This mouthwash is extremely effective in protecting the oral mucosa. You can find out more on page 10. Our products for babies (page 26) also represent a valuable contribution to health protection – starting with newborn babies. As you'd expect, our soothers relax babies. But they can also slow down and even reverse incipient crossbite developments. In addition, we now introduce an oral care product line for young children, with a toothpaste based on enzymes to support the natural, antibacterial effect of saliva and a toothbrush for an effective, smooth and atraumatic brushing experience.

And we have great news for toothbrush fans who love a gentle, effective brushing experience. Now you can try sensationally gentle with the new CS 12460 velvet. 12,460, ultra-fine filaments gently envelop the tooth and sulcus, cleaning with unprecedented care and efficiency. Try it out!

In conclusion dear dentist and dental professional, we support you with products and techniques that are fun to use – and quickly bring visible results. Which encourages your patients to care more for their teeth and gums. This, combined with your incredible efforts to motivate them, means we can both feel proud of our contribution to good health. Thank you for our valuable partnership.

Sincerely,
Ueli Breitschmid

HOW TO ENSURE A BALANCED MICROBIOME: THE LATEST STUDIES

IT STARTS WITH BRUSHING

A young, cavity-free adult's interdental biofilm contains pathogens that can cause approximal cavities. Many of these pathogens cause heart and other systemic diseases. Daily disruption of the interdental biofilm is crucial.

Bourgeois D, David A, Inquimbert C, Tramini P, Molinari N, Carrouel F (2017): Quantification of carious pathogens in the interdental microbiota of young caries-free adults. *PLoS One*, 2017 Oct 10;12(10):e0185804. doi: 10.1371/journal.pone.0185804. eCollection 2017. <https://www.ncbi.nlm.nih.gov/pubmed/29016613>

ORAL CLIMATE CHANGE

As we all know, sugar encourages cavity-causing bacteria. But it also inhibits the growth of many bacteria beneficial to the health of dental enamel... simply because the metabolism of the bad bacteria causes the pH levels in the mouth to decrease.

Peterson S N, Snesrud E, Liu J, Ong A C., Kilian M, Schork N J, Bretz W (2013): The Dental Plaque Microbiome in Health and Disease. Published: March 8, 2013. *PLoS ONE*, <https://doi.org/10.1371/journal.pone.0058487>

PROBIOTIC PASTILLES

If periodontitis patients get into the habit of sucking probiotic pastilles alongside measures such as scaling, the depth of their pockets decreases.

Mitic K, Kaftandzieva A, Popovska M, Ivanovski K, Pandilova M, Georgieva S, Pesevska S, Stojanovska A, Kapusevska B, Janev E, Mijoska A (2017): Probiotics and Oral Health. *Research Journal of Pharmaceutical, Biological and Chemical Sciences*, 8, 1021-9. https://www.researchgate.net/publication/317064037_Probiotics_and_Oral_Health

HEALTHY EATING

How does a low carb diet affect gingivitis and periodontitis? In Woelber et al's study, a group of participants changed to a diet low in carbohydrates but rich in Vitamin C and D, fibre, anti-oxidants and Omega 3 fatty acids. Four weeks later, they had significantly lower inflammation levels.

Woelber JP, Bremer K, Vach K, König D, Hellwig E, Ratka-Krüger P, Al-Ahmad A, Tennert C. (2016): An oral health optimized diet can reduce periodontal and gingival inflammation in humans – A randomized controlled pilot study. *BMC Oral Health*. 2016 Jul 26;17(1):28. doi: 10.1186/s12903-016-0257-1. Erratum in: *BMC Oral Health*. 2016 Oct 6;16(1):109. PMID: 27460471

THE IMPORTANCE OF THE ORAL MICROBIOME

As Kilian et al point out: According to some studies, the bacteria in the mouth closely corresponds with a range of systemic diseases, such as cardiovascular disorders, rheumatic arthritis, strokes, inflammatory diseases of the internal organs, colorectal cancer, respiratory tract inflammation, meningitis and many more. A balanced oral microbiome is highly recommended.

Kilian M, Chapple ILC, Hannig M, Marsh PD, Meuric V, Pedersen AML, Tonetti MS, Wade WG, Zaura E (2016): The oral microbiome – an update for oral healthcare professionals. *British Dental Journal*, 2016, pp. 657-66. DOI: 10.1038/sj.bdj.2016.865.

SHOCKING P GINGIVALIS RESULTS

A study found that 19% percent of participants had porphyromonas gingivalis in their interdental biofilm. These results were shocking – all the participants were young adults with good oral health. Porphyromonas gingivalis can cause loss of jaw bone – and, together with tannerella forsythia and treponema denticola, it can lead to systemic disease

Carrouel F, Viennot S, Santamaria J, Veber P, Bourgeois D (2016): Quantitative Molecular Detection of 19 Major Pathogens in the Interdental Biofilm of Periodontally Healthy Young Adults. *Front Microbiol*. 2016 Jun 2;7:840. doi: 10.3389/fmicb.2016.00840. eCollection 2016

CURAPROX

ORAL HEALTH MEETS COLOR



MICROBIOME 4

PERIO PLUS 8

ORTHO 14

SONIC TOOTHBRUSHES 16

BAD BREATH 24

CURAPROX BABY 26

CURAPROX KIDS 32

PREGNANCY 36

TOOTHBRUSHES 38

IMPLANTS 52

INTERDENTAL BRUSHES 54

PERIODONTITIS 68

TOOTHPASTES 70

WHITENING 78

CURADEN ACADEMY 80

CURADEN CORP. 83



CHLORHEXIDINE WITH A NATURAL BOOST

With 'Perio plus', you get an all-new range of mouthwashes with chlorhexidine and Citrox® natural bioflavonoids extracted from bitter oranges.

The combination of CHX and Citrox® is more powerful than CHX alone. Citrox® has been proven to enhance the efficacy of CHX. With minimal adverse side effects and a great taste and it won't irritate the taste buds.

The best conditions for treatment compliance.

THE NEW POWER COUPLE: CHX AND CITROX®

Chlorhexidine has been prescribed by dentists for decade because it works. But patients often complain that traditional formulations are unpleasant to use. Now you can recommend 'Perio plus'. It tastes amazing and it won't affect the sense of taste. It's also exceptionally effective. The antibacterial qualities of CHX are strengthened by Citrox®, which contains nine natural bioflavonoids extracted from bitter oranges. It's been proven that CHX and Citrox® are more effective than CHX alone.



**EXCELLENT
IN TASTE**

Perio plus

- Treatment compliance done right thanks to a pleasant flavour: developed by Givaudan
- CHX in proven concentrations for mouthwashes, gel and toothpaste
- In combination with Citrox®: the natural bitter orange extract for increased antibacterial effects and the amino acid polylysine for substantivity
- With the copolymer VP/VA for an extra strong protective film on the mucosa and teeth
- With xylitol for a pleasant sweet taste, protection against caries and strong teeth
- Particularly powerful and effective against a very broad range of microorganisms
- Controls and delays the build-up of biofilm (plaque regrowth)
- No impact on the sense of taste
- The discolouration typical of CHX is minimal if used properly
- No alcohol

Antibacterial mouthwashes are sometimes necessary. But to keep this to a minimum, it is important to turn to mechanical cleaning using reliable products. This may sound hard, but it's actually soft and gentle – even on the approximal surfaces; see the 'CPS prime' on page 60.



Chlorhexidine – tried and tested

It has transformed dentistry, especially when it comes to periodontology and oral surgery. It is feared by germs like Porphyromonas gingivalis. For over 50 years, CHX has been the gold standard in dentistry despite its unpleasant taste, enamel discolouration and disruption to the palate. Although these side effects are reversible, it does mean that treatment compliance could be improved.

Bitter orange bioflavonoids

Antiseptic effects dating back millions of years: Citrox® consists of nine bioflavonoids from bitter orange, also known as Citrus aurantium amara. Citrox® antiseptic effect is stronger and broader than that of CHX. It is also antifungal, antiviral, completely natural and free of side effects. In combination with polylysine amino acids, Citrox® retains substantivity like CHX.

CHX with the power of nature: Perio plus

What is more potent, CHX or Citrox®? Malic et al. investigated this question in 2013*. The results are staggering: Citrox® is shown to be more effective than CHX. However: The combination of CHX and Citrox® shows the best performance of all, which is exactly why 'Perio plus' combines CHX with Citrox®. Thanks to this enhanced efficacy, 'Perio plus' can, if appropriate, be dosed at a lower CHX level.

*Malic S, Emanuel C, Lewis MAO and Williams DW. Antimicrobial activity of novel mouthrinses against planktonic cells and biofilms of pathogenic microorganisms. Microbiol Discov. 2013; 1:11. <http://dx.doi.org/10.7243/2052-6180-1-11>

Protective film, anticaries

The copolymer VP/VA has proven to be particularly mucoadhesive and therefore highly effective. The birch sugar xylitol provides both a sweet taste and protection against tooth decay. This is why we use both VP/VA and xylitol in 'Perio plus'.

Prevention and patient loyalty tip

Recommend that patients avoid consuming tobacco and highly discolouring foods during treatment: coffee, red wine, juices, green and black tea, berries, cherries, saffron and curry. Thanks to the strong effect of 'Perio plus', you can, if appropriate, lower the dosage or even reduce the duration of use. Both options further limit discolouration. And, of course, the additional mechanical cleaning also helps minimise discolouration.

What is Citrox®?

Citrox® is a natural bioflavonoid extract of bitter orange (Citrus aurantium amara). These bioflavonoids have a strong – and broad – antibiotic effect, and are antifungal and antiviral.

Over the last few decades, Citrox® has established itself as a preservative in the food industry as well as in surface applications, including for the disinfection of medical facilities.

This has attracted the interest of a number of dental professionals, leading to the 2013 discovery by Malic and her team: Citrox® is even more effective than CHX. There was also a second finding: In combination with CHX, Citrox® has an even stronger effect. Conclusion: Citrox® has been proven to increase the effect of CHX.

Citrox® is a real winner. And it's completely natural.



'PERIO PLUS' IN TRIED AND TESTED CHX CONCENTRATIONS

'Perio plus' is the new CHX range from Curaprox. It combines Citrox® with CHX in the tried and tested concentrations of 0.05 %-0.5 % chlorhexidine gluconate. Introducing 'forte', 'protect', 'balance' and 'regenerate' mouthwashes, 'focus' gel and 'support' toothpaste. Citrox® and CHX: Fighting harmful oral bacteria together.

With hyaluronic acid: 'Perio plus regenerate', 'Perio plus focus gel' and 'Perio plus support toothpaste'. Hyaluronic acid is especially beneficial for the healing of wounds: hyaluronic acid moisturises the upper layers of tissue, promoting growth and stabilisation of the tissue, reducing the risk of complications and accelerating the healing process. This is because hyaluronic acid facilitates cell migration and proliferation, promotes the formation of granulation tissue, combats inflammatory processes and accelerates reepithelialization.

'Perio plus regenerate' mouthwash with hyaluronic acid can be used for up to one month, thanks to only 0.09% CHX. 'Perio plus focus' gel with 0.5% CHX - which also contains Citrox® and hyaluronic acid - is suitable for precise local applications for up to seven days.



Perio plus forte mouthwash

Active ingredients

- Chlorhexidine 0.2%
- Citrox®
- Xylitol
- VP/VA
- Polylysine

'Perio plus forte' offers the strongest effect on plaque regrowth, e.g. for short-term and highly effective intensive care before and after invasive dental treatment.

Perio plus protect mouthwash

Active ingredients

- Chlorhexidine 0.12%
- Citrox®
- Xylitol
- VP/VA
- Polylysine

'Perio plus protect' moderate dosage provides phased antiseptic support for oral hygiene, which is ideal during the treatment of gingivitis, periodontitis and peri-implantitis.

Perio plus regenerate mouthwash

Active ingredients

- Chlorhexidine 0.09%
- Citrox®
- Hyaluronic acid
- Xylitol
- VP/VA
- Cyclodextrin
- Polylysine

Thanks to the addition of hyaluronic acid, 'Perio plus regenerate' accelerates tissue renewal, particularly during the healing process.

Perio plus balance mouthwash

Active ingredients

- Chlorhexidine 0.05%
- Citrox®
- Sodium fluoride
- Xylitol
- VP/VA
- Polylysine

'Perio plus balance' reduces the risk of caries and gingivitis for wearers of orthodontic appliances; in cases of treatment-related dry mouth; in cases of motor impairment. For long-term use.

Perio plus focus gel

Active ingredients

- Chlorhexidine 0.5%
- Citrox®
- Hyaluronic acid
- Xylitol
- VP/VA
- Cyclodextrin
- Polylysine

With the highest dosage of CHX, 'Perio plus focus' inhibits inflammation caused by implants and dentures, and protects against periodontitis and peri-implantitis. It also contains hyaluronic acid to ensure accelerated tissue renewal.

Perio plus support toothpaste

Active ingredients

- Chlorhexidine 0.09%
- Citrox®
- Sodium fluoride
- Hyaluronic acid
- Xylitol
- VP/VA
- Polylysine

'Perio plus support' is free of sodium laureth sulphate (SLS), thus ensuring the success of CHX treatments. At the same time, it contributes to the treatment's success with its own active ingredients.

ORTHO



FOR PATIENTS WITH BRACES

Removing plaque becomes more difficult for patients who wear braces. At Curaprox, we have everything your patients need to give their braces and their oral health the ultimate in tender, loving care.



CS 708 implant | ortho

See page 45
The super compact brush head gives you the perfect cleaning angle. And you can bend the CS 708 to find the exact angle to reach the previously unreachable places.

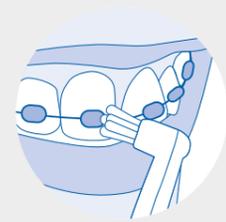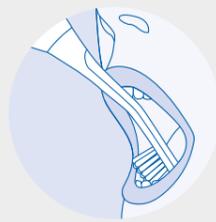

Hydrosonic pro

See page 16
Professional dental care: With 'Hydrosonic pro' you get that professionally cleaned feeling every time you brush - even with braces, implants, sensitive gums... and those hard-to-reach spots.



CS 5460 ortho

- For clean braces, teeth and gums
- So gentle
- 5,460 Curen®-filaments...
- ... each 0.1 mm in diameter
- Special groove in the filaments wraps around braces and clasps
- With practice logo on request
- Developed at the request of dentists



CPS ortho

- For cleaning under wires
- CPS 14 for internal wires
- CPS 18 for external wires
- Our click system for all holders is kinder on the planet

Ortho wax

- Protects oral mucosa
- No more cuts
- Easy to shape
- Easy to apply
- Easy to remove



Perio plus

See page 10
The mouthwash with CHX - and Citrox® natural bioflavonoids extracted from bitter oranges. Reduces the risk of caries and gingivitis for people wearing orthodontic appliances. Also suitable for long-term use.

SONIC TECHNOLOGY INCREDIBLY GENTLE SUPER EFFECTIVE

The **Hydrosonic pro**. Professional teeth cleaning with sonic technology and innovative Swiss brush heads – for exceptional oral care.

Three brush heads developed by oral health professionals: **Power, sensitive and single**. Ultrafine Curen®-filaments. The Curacurve® bend. Imagine enjoying that professionally cleaned feeling every time you brush – even with braces, implants, sensitive gums... and those hard-to-reach spots.

And the **'Black is white' sonic toothbrush** with its 'Carbon' brushhead – featuring activated carbon filaments to further enhance the whitening effect of the 'Black is white' toothpaste.

HYDROSONIC PRO

INCREDIBLY GENTLE, BETTER RESULTS

We developed the 'Hydrosonic pro' sonic toothbrush to give your patients professional support for their oral health. This includes

1. tender loving care for the sulcus
2. preventative and special care for patients with braces or implants
3. the treatment and prevention of gingivitis and periodontitis.

And it's so simple to use – just one button for all functions... with blue LEDs to indicate the function. In seven steps it switches from 44,000 to 84,000 brush movements a minute – to give patients a wonderfully gentle, wonderfully powerful clean. Innovative Swiss brush head technology with Curen®-filaments and Curacurve® ergonomics make the Hydrosonic pro unique. Available exclusively in dental surgeries.



Hydrosonic pro

- Professional teeth cleaning. Every day.
- Tender, loving care for the teeth, gums and sulcus
- Also perfect for patients with braces or implants
- Treats and prevents caries, gingivitis, periodontitis and peri-implantitis
- From 44,000 to 84,000 wonderfully gentle brush movements a minute
- Seven different cleaning levels
- Supports the hydrodynamic effect
- Three innovative brush heads: power, sensitive and single
- 60 minutes brushing from one charge
- The power remains exactly the same throughout the charge
- Charging via USB

Brushing braces and implants with the seven-tuft brush? See page 45.

Larger cleaning area:
improved performance even in
those hard-to-reach spots

Envelop each individual
tooth – innovative water
drop-shaped design



Ultra-fine Curen®-filaments. Gentle.
Effective.

No unpleasant vibrations – thanks to
the rubberised back



Cleans even in hard-to-reach areas... thanks to Curacurve® –
small bend, big difference

Help! Oral health – the fab four

The 'Hydrosonic pro' helps all your patients. It strengthens oral health both with everyday prevention and where careful care is required.

The sulcus: needs gentle, effective care to prevent gingivitis. Ultra-fine filaments and rapid movement reliably remove bacteria from the sulcus and clean without injury.

Implants: The compact brush heads clean hard-to-reach places around implants and fixed prostheses. The single brush head – for tooth-by-tooth and gumline brushing – gives your patients the best prevention against peri-implantitis.

Ortho: The 'Hydrosonic pro' reaches the critical areas with the highest risk of plaque – contact points between teeth and braces, brackets or retainers... and cleans gently and efficiently.

Prevention: Your patients can enjoy exceptional prevention as part of their daily oral health routine... whether they have dental problems or not. The full, gentle power of sonic technology and hydrodynamics cleans thoroughly and cares for the gums, teeth and sulcus... for that professionally cleaned feeling.

A BRUSH HEAD FOR EVERY SCENARIO

Click on and start brushing – safe in the knowledge you’re doing the best you possibly can for your oral health. This is thanks to three innovative brush heads – developed and manufactured in Switzerland – that work gently and effectively. Thanks to Curen®-filaments, brushing feels incredibly gentle – and the droplet shape of the brush heads adds to the efficiency. And thanks to Curacurve® ergonomics, the brush heads easily get to all those hard-to-reach spots.



Brush head ‘power’

- Intensive-yet-caring tooth cleaning
- For use with healthy gums
- With Curen®-filaments
- Droplet-shaped design
- Extra small, extra precise
- With Curacurve® – small bend, big difference
- Reaches all critical areas
- Rubberised back
- Supports the hydrodynamic effect
- Developed with Professor Dr Ulrich P. Saxer
- Made in Switzerland



Brush head ‘sensitive’

- Very gentle on sensitive gums
- For sensitive gums, gingivitis, periodontitis and peri-implantitis
- With Curen® filaments
- Droplet-shaped design
- Extra small, extra precise
- With Curacurve® – small bend, big difference
- Reaches all critical areas
- Rubberised back
- Supports the hydrodynamic effect
- Developed with Professor Dr Ulrich P. Saxer
- Made in Switzerland



Brush head ‘single’

- Extra small brush head
- For the gaps between teeth and braces, brackets, arches, wires, bridges and implants
- For intensive sulcus care
- For cleaning around malpositioned teeth
- For polishing tooth surfaces
- With Curen® filaments
- Extra small for highly precise cleaning
- With Curacurve® – small bend, big difference
- Reaches all critical areas
- Rubberised back
- Developed with oral health expert Jiří Sedelmayer
- Made in Switzerland



The filaments

Compared with standard nylon bristles, Curen®-filaments have the benefit of not adsorbing water. Performance is constant. Curen® stays incredibly soft and feels wonderful – even at 84,000 pulsations per minute.



The ergonomics

Curacurve® is a bend in the brush of only 10 or 15 degrees. But it makes a world of difference. Now you can easily clean the difficult areas – the back of the front teeth, behind the molars and premolars.



Enveloped
Thanks to Curacurve® technology, the filaments reach every spot and clean perfectly. Prof. Ulrich P. Saxer designed the brush heads ‘sensitive’ and ‘power’ in the shape of a water droplet.

They feature long filaments on the outside and short filaments on the inside. This allows each tooth to be enveloped individually for optimum mechanical and hydrodynamic cleaning.

Comfort

Electric toothbrushes can feel unpleasant if the back of the brush head touches a tooth. So we rubberised the back of the brush head... in a way that doesn’t impair cleaning performance. That’s innovative Swiss brush head technology.



THE WOW FACTOR! DENTISTRY MEETS DESIGN

The usual gentle, effective, cleaning power of Curaprox. Super styled. Treat your teeth to up to 84,000 wonderfully gentle brush movements a minute. The Hydrosonic 'Black is white' even cleans the interdental spaces thanks to the hydrodynamic effect. Also perfect for patients with implants, bridges and braces.



Hydrosonic 'Black is white'

- The wow factor! Dentistry meets design
- 44,000, 64,000 or 84,000 wonderfully gentle brush movements a minute
- Brush heads with activated carbon filaments for gentle whitening
- More than 60 minutes of brushing
- Charging via USB
- The hydrodynamic effect
- Developed in Switzerland with Prof. Ulrich P. Saxer

The hydrodynamic effect

Now your patients can enjoy the powerful results of water in motion. The brush head swirls around the mixture of saliva, water and toothpaste at high speed. There's no hiding place for the bacteria and plaque. It's instantly flushed out. That's the hydrodynamic effect. To see the difference... Hold the brush horizontally in a glass of water, 2mm-3mm under the surface. Watch the water bubble. Then repeat the test with a standard electric toothbrush.

CURAPROX
BLACK IS WHITE





Curen®-filaments with activated carbon particles

Droplet-shaped arrangement

Rubberised back




Curacurve® ergonomics

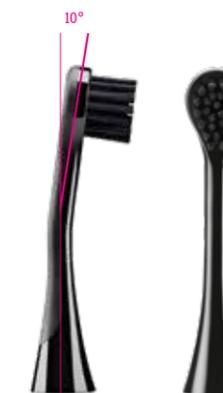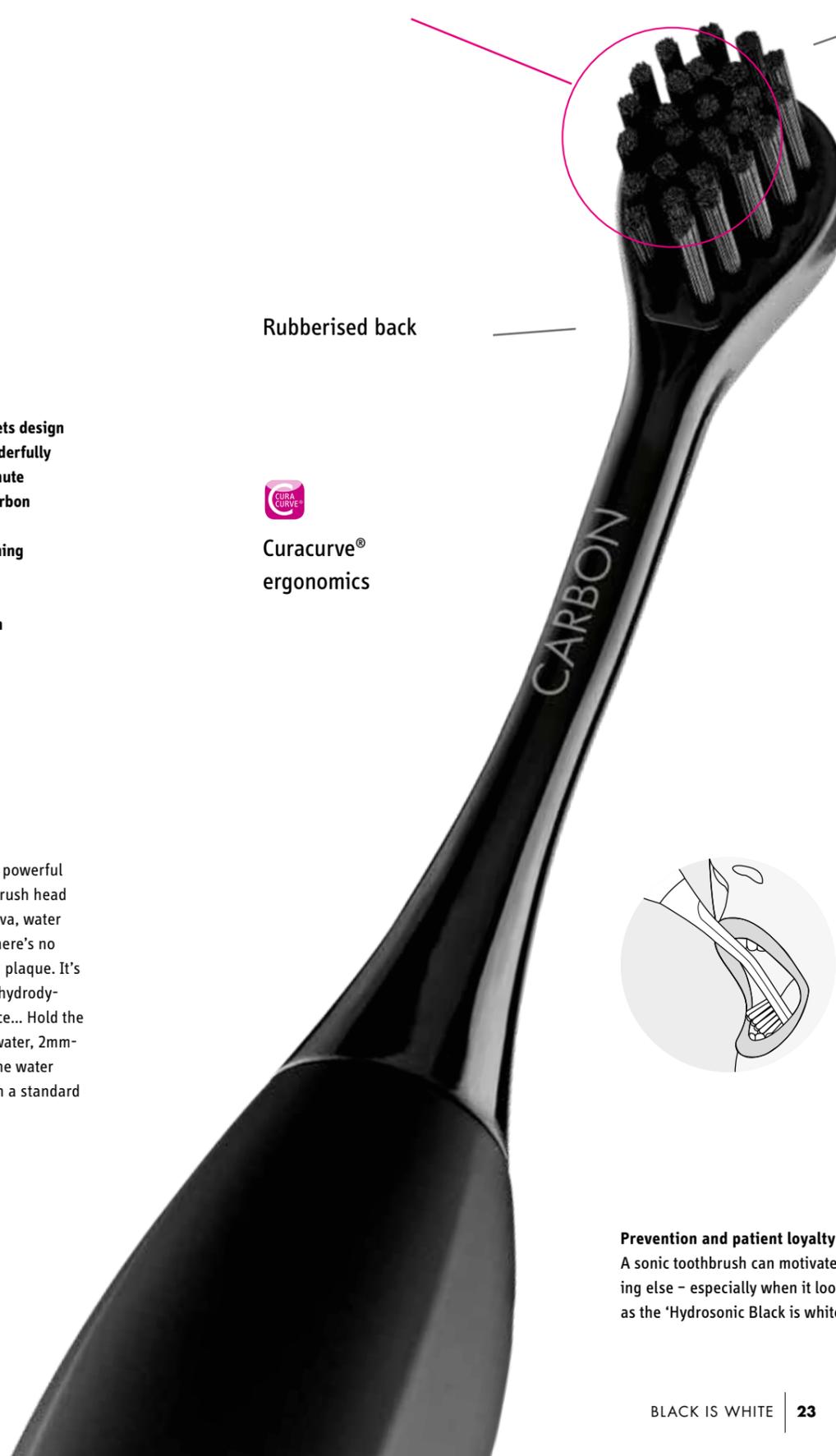

Anything else?

There's the Curacurve® ergonomics: this bend in the brush head makes all areas in the mouth easily accessible. And with its black and stylish finish, it perfectly matches the Hydrosonic 'Black is white'.

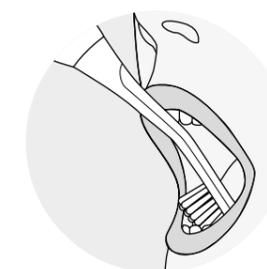

Prevention and patient loyalty tip

A sonic toothbrush can motivate like nothing else – especially when it looks as good as the 'Hydrosonic Black is white'.

'Black is white' is the new black
'Black is white' products have innovative technology at their core. And that's where the technology should stay – inside the product. On the outside, oral health products should look and feel beautiful.

Your eyes help you brush...

... and that's why we're obsessed with design as well as efficiency. Design matters. It motivates your patients to brush regularly. And it gives them stylish colourful or black-and-white accents in the bathroom.

Form and function

But we will never compromise on functionality and performance. This is how we developed our 'Black is white' concept – form and function in perfect harmony.

BAD BREATH

FOR PATIENTS WITH BAD BREATH

When it comes to feeling confident about teeth, gums and oral health, there are two main factors to consider – especially for younger people. The first is feeling great about an attractive smile. The second is feeling self-assured about breath. For younger people, the tongue is often the main cause of halitosis. For older people, periodontitis and reduced saliva production can also be major factors. Careful oral hygiene is the best way to avoid or eliminate halitosis – and to feel fresh and confident. Halitosis affects 25% of the population. So your recommendations to your patients can make a significant difference.



Hydrosonic pro

See page 18

The bacteria and plaque that cause bad breath are instantly flushed out as the brush head swirls around the mixture of saliva, water and toothpaste at high speed.



Enzycal

See page 72

Gentle, mild, effective, fresh-tasting toothpaste. With enzymes that boost the protective power of saliva... helping your patients to limit the growth of the bacteria that cause bad breath – and to prevent tooth decay.



CTC 202 tongue scraper

Features two blades – one smooth, one serrated – which work together to thoroughly remove the bacteria that build up on the tongue. Fresh breath confidence in two seconds!



CPS prime

See page 60

Super-fine filaments for super-clean interdental spaces. The foundation for great oral health and fresh breath. Help your patients to say no to gum inflammation and bad breath. It's easy. Once in and out. Done.



Perio plus

See page 10

For longer-term cases of bad breath, try 'Perio plus balance' – a naturally effective mouthwash. It combats bacteria by combining CHX and Citrox® natural bioflavonoids extracted from bitter oranges.



'Black is white' chewing gum

Temporary relief from bad breath with a revolutionary enzymatic whitening chewing gum... that actually starts a flow of saliva. Chew your way to a healthier mouth and whiter teeth.



ORAL HEALTH FROM DAY ONE

Our ultimate **soother** supports – instead of hinders – the normal development of gums, teeth, the palate and jaw... and ensures optimal breathing.

The teething ring – with integrated training toothbrush – helps your child prepare for the day they start brushing their teeth.

And the **Curaprox baby toothbrush** gives the gentlest care to precious teeth.

CALMS BABIES AND THEIR PARENTS

The Curaprox baby soother with its flat suction tip, side wings and silicone membrane looks different from all the others. Developed by orthodontist Dr Herbert Pick, this is the ultimate orthodontic soother. It relaxes babies without the risk of dental displacement or the danger of jaw deformation. It prevents the teeth from moving and prevents difficulties with overall health. Our soother supports the natural development of a baby's teeth, gums, palate and jaw... and supports language development and the natural swallowing process. And it ensures optimal breathing through the nose – good for concentration and restful sleeping. And, as such, it calms babies – and their parents!



Biofunctional soother

- The ultimate orthodontic dummy
- Supports the natural development of your baby's teeth, gums, palate and jaw...
- ... and ensures optimal breathing
- Natural breast sensation supports the natural swallowing process
- Developed by distinguished Dr Herbert Pick
- Biofunctional
- No harmful substances
- 3-14 kilograms
- 4 colours, 3 sizes
- With steriliser and transport box
- Recommended accessory: Soother holder

The flat tip helps to ensure the palate doesn't take on an abnormal, high-arched shape. This prevents teeth from moving and reduces problems with breathing, sleep and concentration.



The side wings guide the suction pressure to the jaw rather than the palate, ensuring the dental arch develops a nice, round shape... giving the milk teeth – which act as the placeholders for adult teeth – the space they need.



The special, soft silicone membrane creates the same, pleasant sensation as a natural nipple and supports the natural swallowing process.



Just click!
The soother holder stops the soother falling on the floor – where it can get dirty and contaminated. It's simple and secure. It won't get lost. It's always there for your child.

Free from harmful substances: No Bisphenol A (BPA), no Phthalate, no azo-dyes.



Tip
The cleaning and transport box: Practical. Essential.



Top tip

Give away a free Curaprox baby soother to your pregnant patients. It's a wonderful gesture that your patient will remember for a long time – and it encourages the best oral health from day one.

Are soothers good for your child?

When teenagers need orthodontic treatment, it's often because they used a dummy – or sucked their thumb – as a young child. 'In western Europe, 70-80% of children have a dummy and one in two of them develop a crossbite or an open bite,' says Dr of Dental Science Herbert Pick. In fact, only 20% of such 'malocclusions' are genetic.

Children with such problems are more prone to mouth breathing... which, in turn, leads to an insufficient supply of oxygen to the brain. This can cause developmental disorders. Dr Pick, who developed the Curaprox soother, says: 'There is a clear link between jaw abnormalities and a tendency towards mouth breathing, which can lead to problems such as disturbed sleep and difficulty concentrating. Poor posture has also been cited in the literature as a consequence of mouth breathing, along with speech and concentration disorders and even hyperactivity and heart conditions.'

For Dr Pick, a distinguished orthodontist in Kufstein, Austria, developing the ultimate soother became a personal mission... after his son developed a crossbite. 'I felt challenged, not just as a father but also as an orthodontist, and I wanted to know exactly why this crossbite happens and how it can be prevented,' he says.

As with all our baby products, the Curaprox soother was developed using our 'biofunctionality' ethos... thinking about cause and effect – and creating products that promote growth in a natural way.



TAKE AWAY TEETHING PAIN

For babies, the teething ring is a toy. For parents it takes away teething pain – their baby plays happily... and they're content. Developed by Dr of Dental Science Nadja-Marina Kellerhoff, the Curaprox Baby teething ring promotes the development of sensory and motor skills. It also gets babies used to brushing in a playful way thanks to the integrated massage and teaching brush. Teething is painful. But there's no need for parents to reach for the ointment. Babies will instinctively bite down on the teething ring and it will take their pain away. What a relief!



- **Teething ring top tip 1**
When the first tooth breaks through, it's a celebration for parents... but not for the babies. They're in pain. Tell parents to try the Curaprox Baby teething ring before they go reaching for pain relief gel.
- **Teething ring top tip 2**
Keep the teething ring in the fridge to intensify the child's relief!

Teething ring

- Makes teething more pleasant – for babies and their parents!
- Promotes the development of sensory and motor skills
- A rattle to awaken the senses
- With massage and teaching brush
- Biofunctional
- Developed by Dr Nadja Marina Kellerhoff
- 0-24 months
- 3 colours
- No harmful substances



Oral health from the very beginning – thanks to biofunctionality

Parents naturally want the best for their children. And we want the best for all our customers young and old. All Curaprox Baby products are developed with our biofunctionality ethos. That means analysing cause and effect – and designing products to support growth naturally. A biofunctional approach crucially allows us to recognise long-term consequences such as breathing problems. Our products are designed to eliminate undesirable effects. Put simply, biofunctionality means our products support – instead of hinder – the natural development of babies and small children.

TENDER LOVING CARE

The Curaprox Baby toothbrush, developed by Dental Professor Adrian Lussi, makes cleaning teeth fun and joyful and encourages quality parent-child time. The compact brush head reaches all areas of the mouth – ensuring every tooth is cleaned properly. At the same time, the rubberised brush head protects the oral mucosa. The 4,260 densely-packed Curen®-filaments – each with a diameter of just 0.09 mm – give gentle, exceptional cleaning. Brushing becomes more like stroking. The rounded handle ensures children get into the habit of holding their toothbrush correctly – with gentle pressure. And it fits perfectly into the hands of both children and parents. It goes without saying that it's toxin-free and completely safe to use.

Top tip for patients
Try two baby toothbrushes: One for the child to take ownership of and to playfully practice brushing with. The other for parents to use afterwards to ensure a thorough clean. Find more tips and tricks at curaprox.com/baby

Baby toothbrush

- Tender loving care
- Biofunctional
- Developed by Professor Dr Adrian Lussi
- 4,260 gentle Curen®-filaments with a diameter of 0.09 mm
- Rubberised brush head
- Round handle
- 3 colours
- Free from harmful substances





«HAVE YOU BRUSHED YOUR TEETH YET?»


Trying to motivate the little ones to brush their teeth can become exhausting – for them and for their parents. Let's say goodbye to endless nagging.

But how?

Firstly, the toothpaste tastes so good. The strawberry, watermelon and mint flavours are a taste sensation.

Secondly, the toothbrush feels so pleasant. It's gentle. And it looks bright and beautiful.

Parents may ask about the best toothbrushes and toothpastes for their little darlings. Reassure them that Curaprox products are designed in collaboration with dental professionals.

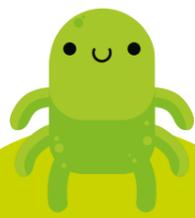

No matter whether a child's permanent dentition is partially or fully erupted, our kids toothbrush ensures that the surrounding gum is uninjured during brushing, enabling children to effectively remove plaque that causes tooth decay and gingival bleeding.

Prof. Dr. Isabel Martínez Lizán



All those qualities embodied by the kids toothbrush provide an effective, smooth and atraumatic brushing experience for teeth and gums – leading to a high acceptability among children.

Dr. Med. Dent. Franka Baranović Huber

GENTLE CARE FOR PRECIOUS TEETH



Children's toothbrush

- Incredibly gentle thanks to 5,500 filaments – each 0.09mm in diameter
- Filament length optimised for children at 7.5mm
- Extra small, compact, slightly angled brush head
- Extra slim octagonal handle – perfect for any brushing technique
- Comes with a protective cap
- Recommended for children aged 4-12
- Six joyful colours
- Made in Switzerland

Father Christmas. The Easter Bunny. The Tooth Fairy. Eventually, kids may start to doubt the magic.

But one thing is for certain: The perfect toothbrush for children exists. Help children to begin the lifelong habit of cleaning gently and thoroughly with a brush that features a multitude of ultra-fine filaments. Keep decay at bay. And give their oral health the best possible start.

Top tip for patients

The Curaprox kids toothbrush is children's new best friend when they lose their baby teeth and their adult teeth start to come through. Gentle. Effective.

GREAT TASTE. HEALTHY MOUTH.

Finally, a toothpaste that makes kids excited about brushing their teeth. Funky fresh strawberry, watermelon and mint give them a party in their mouth. Enzymes support the natural, antibacterial effect of saliva. As does fluoride – with three options: 950 ppm, 1450 ppm for older children and fluoride free. Now parents can feel secure they're doing everything they can to nurture their child's oral health. No harmful substances. Obviously.



ENZYMATIC PROTECTION

The enzymes amyloglucosidase and glucose oxidase in Curaprox kids toothpaste boost the protective power of your child's saliva to fight harmful bacteria. Better protection. Naturally.

Children's toothpaste

- Great taste. Healthy mouth.
- Strawberry, watermelon and mint flavours
- No potentially harmful SLS
- Amyloglucosidase and glucose oxidase support the natural, antibacterial effect of saliva
- Xylitol has an anticariogenic effect
- Low RDA – protects their tooth enamel
- Two fluoride strengths – 950 ppm from the age of 2 and 1450 ppm from the age of 6
- Also available fluoride free
- Easy-to-squeeze tube. Simple to open and close.
- No ingredients derived from animals
- Made in Switzerland



From the first tooth: Strawberry: 0 ppm fluoride.

From the age of 2: Strawberry: 950 ppm fluoride.

From the age of 6: Mint: 1450 ppm fluoride.

FOR PREGNANT PATIENTS

Pregnant patients can find it more difficult to look after their oral health. Hormonal changes make them more susceptible to a build-up of plaque – and also loosen the connective tissue. This increases the blood supply to the gums... leading to softer, more sensitive gums that swell more easily. That's why it's so important to recommend very gentle aids... which can, nevertheless, effectively remove and disorganise bacteria and biofilm – on the enamel, in the sulcus and in the approximal region.

CPS prime

See page 60

'CPS prime' interdental brushes give exceptional cleanliness in those difficult-to-reach narrow spaces, thanks to super-fine, extra-long, ultra-resilient filaments. Failing to clean the spaces between the teeth is a major health risk for anyone. And studies have also suggested a possible link between gum disease in pregnancy and a higher risk of premature birth. 'CPS prime' help to deal with or prevent gingivitis and therefore the further risk of periodontitis.



Enzycal 1450

See page 72

Enzycal 1450 toothpaste with fluoride helps your pregnant patients to enjoy the best oral health – delivered in the most gentle, mild and effective way. Your patients are more at risk of tooth decay during pregnancy. Sodium fluoride strengthens the tooth enamel and, after brushing, makes it resistant to acid attacks. It also remineralises the enamel, including areas that have come under attack... and minimises the growth of the harmful bacteria that can lead to tooth decay.



Pregnancy cravings & morning sickness

Pregnancy cravings

For many pregnant women, cravings start in the first trimester – as early as five weeks into pregnancy – and get stronger in the second trimester. Those who crave acidic or sweet food or drink should take special care. The acidity of things like cola, orange juice, fruit and sugary sweets attack the structure of tooth enamel. Advise your pregnant patients not to eat or drink too many of these treats! When they do, they should take a sip of water straight afterwards and wait at least half an hour before brushing – the teeth enamel becomes weaker and more at risk of damage straight after consuming acidic food and drink.


Morning sickness

Similarly, with the nausea and vomiting common in pregnancy – often known as morning sickness – the acid in vomit can damage the teeth enamel. Patients should rinse the mouth with water after being sick and wait at least half an hour before brushing.

CS 5460

See page 42

Because gums are soft and even more sensitive during pregnancy, it's important to clean them gently. That's why our iconic CS 5460 toothbrush features 5,460 ultra-fine, gentle Curen®-filaments, densely packed on the compact brush head. This gives the patient gentle but exceptional power to clean the teeth and the gumline without causing any injury.

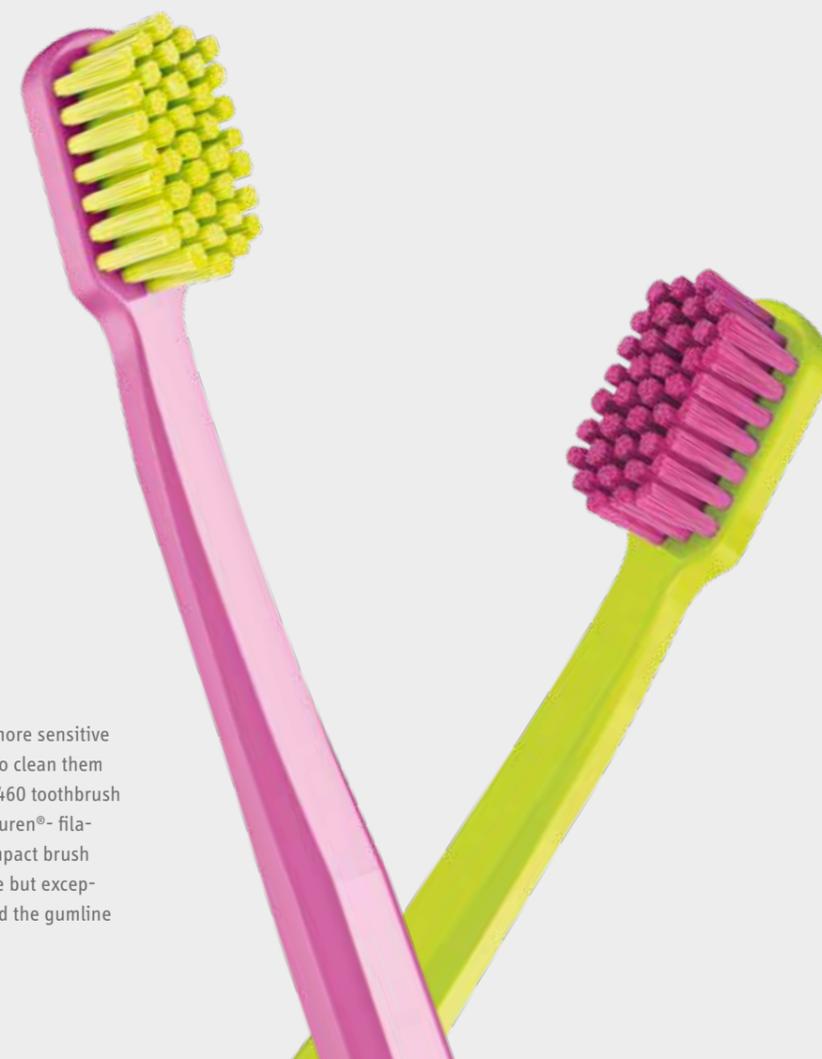



EASE THE PRESSURE

Why do dental professionals recommend our iconic **CS 5460**? Two words. Gentle. And efficient. The precise **CS smart** is even more gentle. And our brand new **CS 12460 velvet** takes things to the next level... from the gentle caress when you start brushing to the silky smooth sensation after cleaning.


Meanwhile, the **CS 1006 single** – with its single tuft – is the ultimate in tooth-by-tooth and gumline cleaning.

Also introducing the fully flexible **CS 708 implant | ortho** with its super compact seven-tuft head to make sure your ‘hard-to-reach’ spots are now simple to get to.

After surgery and during radiotherapy of the oral cavity, the **CS surgical** comes to the rescue with an especially gentle clean.

And you can start them young with the **kids toothbrush** – recommended for children aged 4+.

LOVE YOUR GUMS

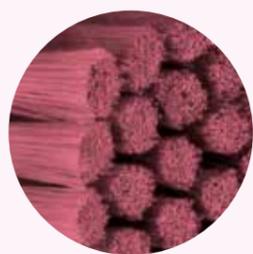

These are incredible 12,460 Curen® filaments. Each filament has a diameter of hair-thin 0.08 mm.

Anyone who cleans with our toothbrushes knows how gentle they are. We've gone one better now: May we introduce you to the new 'CS 12460 velvet'? She comes with a total of 12,460 filaments standing close to close, forming a highly efficient yet caressingly soft cleaning surface. She leaves nothing but cleanliness, strokes sulcus and every single tooth; the enamel remains undamaged. Well, there's just one thing you have to do: Try it out!



NEW

CS 12460 velvet


- The gold standard
- 12,460 Curen®-filaments
- Each 0.08 mm in diameter
- Slightly angled, compact brush head for optimal cleaning
- Octagonal handle to ensure your patients brush at the correct 45 degree cleaning angle
- 36 colour combinations
- Made in Switzerland

36 colour combinations. Matt. Elegant. Beautiful.

Filaments



Brush

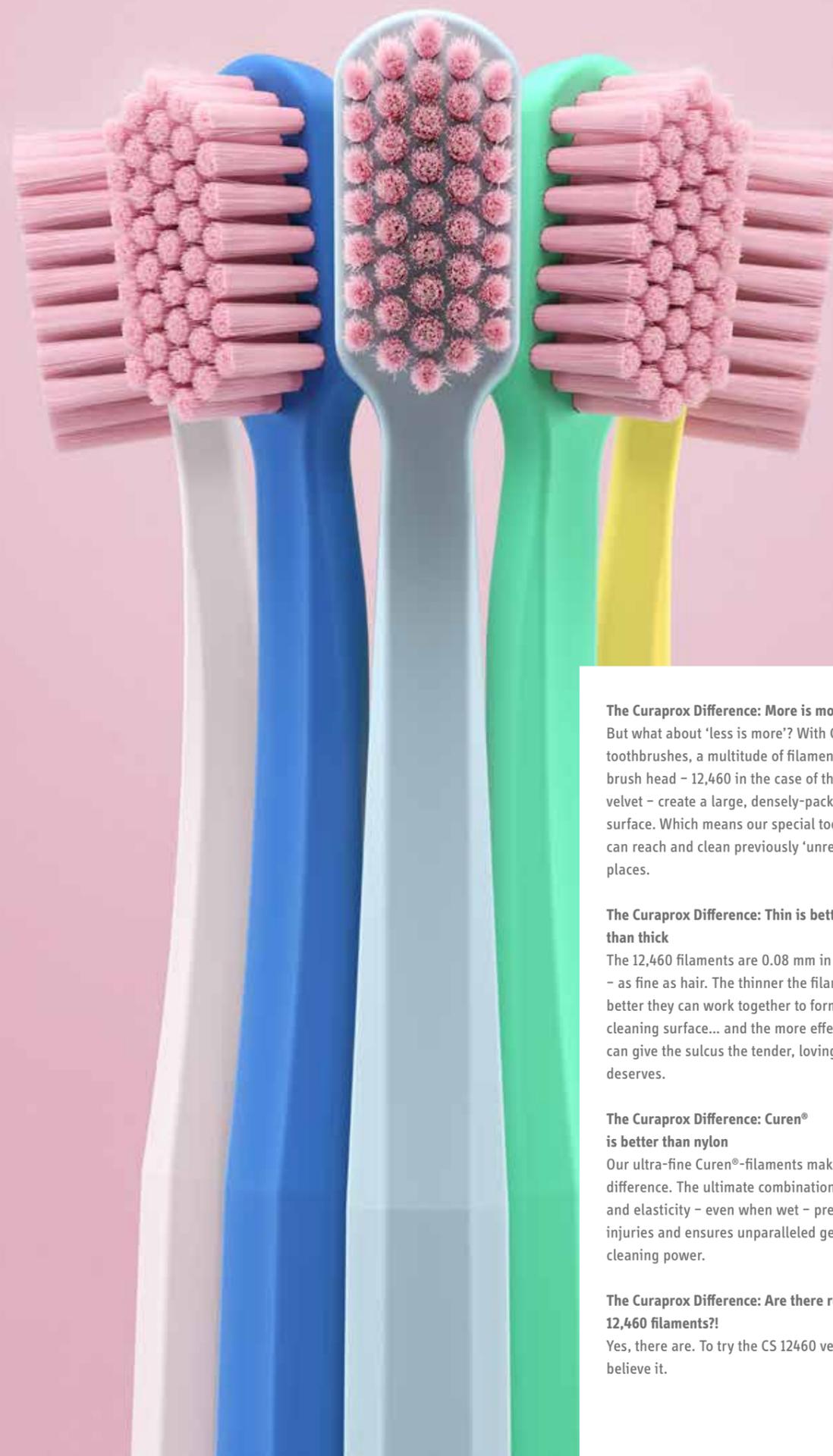

The Curaprox Difference: More is more

But what about 'less is more'? With Curaprox toothbrushes, a multitude of filaments on the brush head – 12,460 in the case of the CS 12460 velvet – create a large, densely-packed cleaning surface. Which means our special toothbrushes can reach and clean previously 'unreachable' places.

The Curaprox Difference: Thin is better than thick

The 12,460 filaments are 0.08 mm in diameter – as fine as hair. The thinner the filaments, the better they can work together to form a dense cleaning surface... and the more effectively they can give the sulcus the tender, loving care it deserves.

The Curaprox Difference: Curen® is better than nylon

Our ultra-fine Curen®-filaments make all the difference. The ultimate combination of fineness and elasticity – even when wet – prevents injuries and ensures unparalleled gentle, cleaning power.

The Curaprox Difference: Are there really 12,460 filaments?!

Yes, there are. To try the CS 12460 velvet is to believe it.

WHEN YOUR GUMS FALL IN LOVE

The CS 5460. Relentless on plaque. Relentlessly gentle on your teeth and gums...

All thanks to 5,460 ultra-fine, Curen®-filaments densely packed on the brush head to give gentle, exceptional cleaning power... so you can remove plaque effectively without causing damage. And what about that iconic design? An explosion of beautiful, joyful colour.



CS 5460

- Brushing your teeth is pure joy
- 5,460 Curen®-filaments...
- ... each 0.1 mm in diameter
- Compact, angled head
- Octagonal handle for 45-degree brushing angle
- 36 colour options
- Special / Limited editions
- With practice logo on request
- Made in Switzerland


The personal touch
We can print your practice logo on the CS 5460. Now, every time your patient brushes they'll be thinking of you!




5,460 is the magic number
5,460 filaments – up to ten times the number of bristles used in conventional brushes – tightly packed together on the brush head. This extra dense cleaning surface is designed to prevent tooth decay and gum disease. The CS 5460 is plaque's worst nightmare.

A multitude of ultra-fine filaments


You see it every day. Abrasions, shrinking gums, bleeding gums, exposed tooth necks. All caused by brushing with hard bristles and too much pressure.

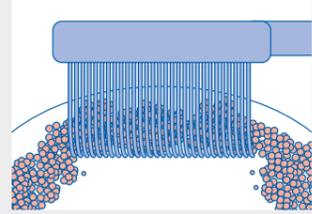

Curaprox toothbrushes are gentle and effective – and remove plaque thoroughly while preventing enamel damage. How do we do this? Curen®-filaments are more rigid than conventional nylon bristles so they remain stable. They don't adsorb water so they stay dry in the mouth. All of which allows us to produce toothbrushes with an incredible amount of super-fine filaments.

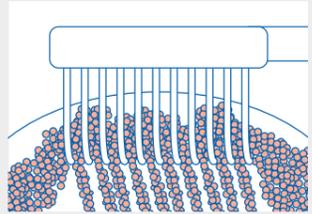

By comparison, a standard toothbrush uses way less, loosely-arranged bristles – leaving lots of plaque on the tooth. The bristles are often so hard that the only way to try to clean in hard-to-reach spots is by using excess pressure. Damage is inevitable.

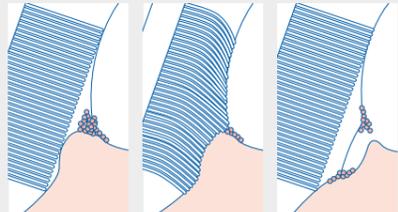

Curaprox toothbrushes all feature an extremely high number of ultra-fine filaments. The CS 5460 boasts 5,460 filaments – each with a diameter of just 0.1mm. The CS smart has 7,600 filaments at 0.08 mm in diameter. And the CS 12460 velvet has 12,460 filaments at 0.08 mm. The result? An extremely dense area of bristles that disrupts and removes plaque considerably more thoroughly and gently.

CS SMART

OH YES!

The CS smart toothbrush – based on the CS 5460 – features an extra small brush head... but more Curen®-filaments. Dental experts recommend the CS smart because it's such a precise brush. Some 7,600 filaments are tightly packed together to give gentle, exceptional, cleaning power. The huge number of filaments make brushing unbelievably gentle and injuries all but impossible.



CS smart

- 7,600 Curen®-filaments
- Each 0.08 mm in diameter
- Extra small brush head
- Octagonal handle
- 18 colour options
- Developed by oral health expert Jiří Sedelmayer
- Made in Switzerland

CS 1006

PASSION

The CS 1006 – with its single tuft – doesn't look like a normal toothbrush. 'Is it even a toothbrush,' you might ask. You'd better believe it. Tooth brushing doesn't get more precise than this – it's the ultimate in tooth-by-tooth and gumline cleaning. Toothbrush aficionados love making those precise, gentle movements along the gumline with the CS single – cleaning in places where it hardly seems possible. Even around implants... and between teeth and braces, brackets or retainers.



CS 1006

- For the teeth cleaning connoisseur
- Single head for sulcus care
- Curen®-filaments – each 0.10 mm in diameter
- Edge length 6mm
- 24 colour options
- Developed specifically for Jiří Sedelmayer's 'solo technique' brushing method
- Developed in Switzerland. Made in Germany.

Want to find out about cleaning your teeth with 44,000 brush movements a minute? See page 18.



The solo technique – the short version


Hold your brush vertically and apply to one side of the tooth's surface. Carefully move the filaments of the brush head so they're just touching the gumline. Now move the CS single with precise, gentle, 'shaking' movements along the gumline – from one side of the tooth to the other.

Go to curaprox.com for full instructions & video guide



PERFECT BRUSHING ANGLE

'I can't get to those hard-to-reach spots in my mouth.' You can now with the CS 708 implant | ortho. Just ask the patients with implants, fixed prostheses, braces and lingual braces who have joined the CS 708 revolution. The super compact brush head gives you the perfect cleaning angle. And you can bend the CS 708 find the exact angle to reach the previously unreachable places.

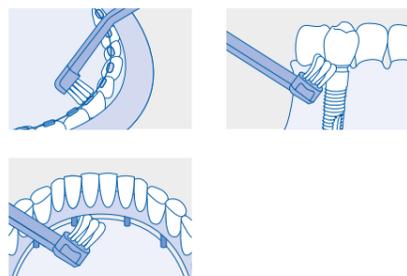

How to use: Simply clean where you need to without pressure: implants, fixed prostheses, wires, braces, lingual braces, brackets. Bend the CS 708 to perfect your brushing angle if required.

Top tip for patients

Make sure your patients understand that, if they brush well, they prevent inflammation... improving overall health.

Thanks to its angled head, this toothbrush is the perfect solution for cleaning especially tight and hard-to-reach places around implants and orthodontic devices.

Many people wear braces to correct dental positioning. But it makes keeping the teeth clean a real challenge. If the patient cleans carefully around and between wires and braces with the right toothbrush, it's possible to keep teeth and gums healthy throughout orthodontic treatment.

The CS 708 has been developed to clean all kinds of oral implants, including fixed prostheses. It's the perfect hygienic tool to help patients avoid mucositis and peri-implantitis.

Professor Dr Kirsten Warrer, University of Aarhus, Denmark, and Professor Dr Hugo Roberto Lewgoy, University of Sao Paulo, Brazil.



CS 708

- The perfect cleaning angle
- Cleans hard-to-reach places
- Perfect for patients with implants, fixed prostheses, braces and lingual braces
- Super compact seven-tuft head
- Curen®-filaments – each 0.12 mm in diameter
- Bendable handle
- 18 colour options
- Developed with Professor Dr Kirsten Warrer and Professor Dr Hugo Roberto Lewgoy
- Made in Switzerland

Now ortho and implant care is easy!
See also pages 14, 64 and 66.

MEDICAL TOOTHBRUSH

Brushing when brushing seems impossible. The CS surgical mega soft is the perfect post-surgery toothbrush. After surgery and during radiotherapy of the oral cavity – especially when the gums are still wounded or particularly sensitive – good oral hygiene is crucial. 12,000 super-flexible, extra-fine filaments come to the rescue. This is extra gentle dental care that guarantees oral hygiene – even during special medical situations. Developed in collaboration with Professor Dr N. P. Lang, Dental Clinics of the University of Berne.



CS surgical mega soft

- Medical toothbrush
- Perfect post-surgery and during radiotherapy
- 12,000 Curen®-filaments
- Each 0.06mm in diameter
- Developed with Professor Dr N. P. Lang
- Developed in Switzerland. Made in Germany.

Support the healing process: with the help of CHX, Citrox®, hyaluronic acid and PVP/VA; see page 10.

Tender, loving care

You can use the CS surgical alone on the operation site for two to five days. This prevents further injury while providing an effective clean. This toothbrush is astonishingly gentle. When you run your finger over the filaments you will notice virtually no resistance at all. Curen®-filaments are the ultimate in tender, loving care.

CURAPROX TRAVEL SET

THE PERFECT TRAVEL COMPANION

Brush your teeth on the move with the ultimate Travel set for oral health care. No compromise on style. No compromise on your oral health. Your favourite Curaprox products in your suitcase, your hand luggage – or in your pocket. This vibrant, stylish travel kit measures just 9.5cm x 6cm.

Wave goodbye to cavities and gum inflammation and say hello to healthy, beautiful teeth!



Curaprox Travel set


- CS 5460 travel toothbrush: Relentless on plaque. Relentlessly gentle on your teeth and gums
- 'Be you' toothpaste: Protects and cares for teeth and gums. Gently and effectively whitens
- 'CPS prime' interdental brushes: Effective brushing in the spaces between your teeth
- Now you can always give yourself the oral health care you deserve – wherever you are in the world.



The Travel set: Your smallest case is now your most important case – no matter where you travel.

SPECIAL EDITIONS

Exclusive artistic interpretations of our iconic CS 5460 toothbrush. Eye-catching bathroom accessories to suit every mood – sold in twin packs. These exquisite designs are the result of collaborations with local, international, upcoming and established artists. Bring infectious, multicolour joy to your oral health routine.

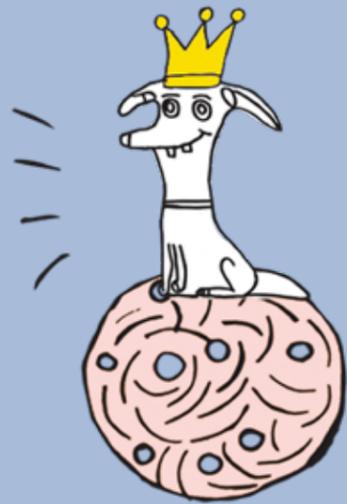

SPECIAL EDITIONS

AVAILABLE WHILE STOCK LASTS



Summer Edition
Cool down with a refreshing ice lolly. Laze on a lilo in the pool. Paddle board in the sea. The toothbrush for summer lovers.



Fun with Curaprox
Smiling while you brush with your stylish, burgundy CS 5460 toothbrush... as you remember uncontrollable laughter at the weekend. Fun times.



Tiger Edition
Courage. Beauty. Power. Strength. Confidence. Style. Everything you need to seize the day. Go get 'em tiger. Rooooooar!



Love Edition
Love is... taking it in turns to brush your teeth... for giving them for leaving the toothpaste open... sharing a toothpaste-fresh kiss.



Little Bacteria Edition
Bright green and bold purple toothbrushes for your kids. Inspired by bacteria to give your little monsters everything they need for brushing fun.

FOR IMPLANTS PATIENTS

You can help your patients to do everything they can to care for their dental implants with the least amount of hassle... reducing the risk of inflammation and peri-implantitis. With the best oral hygiene, your patients' implants can last for decades.

CS 708 implant | ortho

See page 45


The perfect cleaning angle. Cleans hard-to-reach places. Perfect for patients with implants, fixed prostheses, braces and lingual braces.



CPS implant

See page 64

Perfect for cleaning large interdental spaces. Extra-long, extra-fine, extra-rounded filaments. System with five sizes.

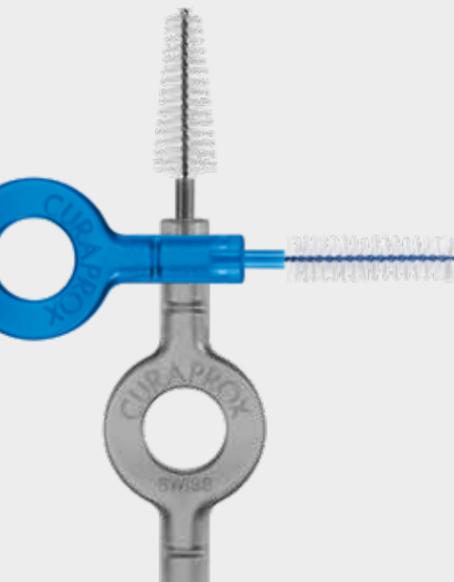

Perio plus

See page 10

Contains CHX together with CitroX® natural bioflavonoids – extracted from bitter oranges – to reinforce the antibacterial effect.

CS 1009

- Like our CS single 1006 – but with a trim that's 3mm longer – it's one for toothbrush aficionados
- Clean millimetre by millimetre, tooth by tooth.
- Optimal care for implants, bridges and braces and isolated teeth – or teeth in an unusual position.
- Curen®-filaments





CS 5460

See page 42

The joy of tooth brushing. Even with implants.


Special care floss

- A dental floss designed for implant care
- Elastic microfibres allow for a novel, efficient cleaning movement around the implant and along the gums
- Tense, relax, super clean
- Gentle, novel, effective implant care

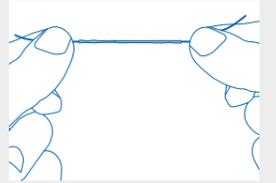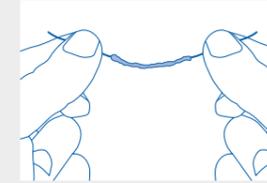

Application

1.

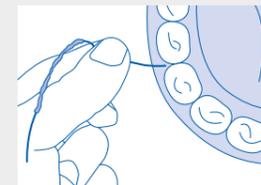


Use the threader to guide between the teeth.

2.



Guide along the sulcus around the implant.

3.

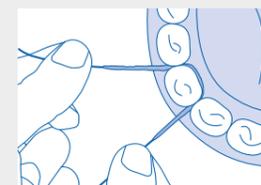

Use two hands to hold each end; clean the sulcus by gently pulling and releasing the implant saver.



THE INTERDENTAL SECRET

Cleanliness in the spaces between the teeth is the best foundation for the oral health of your patients. The **CPS prime** inter dental brush for primary prophylaxis prevents side cavities, gingivitis and periodontitis. The smallest ‘CPS prime’ is just twice as thick as dental floss.

The **CPS perio** for secondary and tertiary prophylaxis was developed for periodontitis patients and is recommended for larger interdental spaces. The **CPS implant**, with its long, fine, rounded filaments is gentle on implants, crowns, bridges and after surgery.

We’re proud of our **CPS system** – measure using **inter dental acces probes** and instruct using the **inter dental acces chart**. Everything’s on hand with the convenient chairside box. Easily motivate your patients by showing them how easy, pleasant and effective it is to use their inter dental brushes.

WORKING WITH THE CHAIRSIDE BOX IN THE DENTAL PRACTICE

The chairside box offers hygienic storage for the Curaprox range of interdental brushes as well as interdental access probes (IAP) and the UHS 415 autoclavable holder.

Your patient is comfortable, so let's begin. First, take an IAP out of the box and attach it to a holder.



1.

Practical and hygienic, directly within your reach.

You can also sell more brushes in your practice. We would be happy to advise you on displays and sales sizes.




7.



Your patient is ready to use their interdental brushes at home: Give them one brush of each size needed. Simply place the brushes inside the flyer that explains how to use them correctly.

Test that the sizes are right. Add the definitive sizes to the interdental access chart – ideally with a pen in the matching colour.

6.

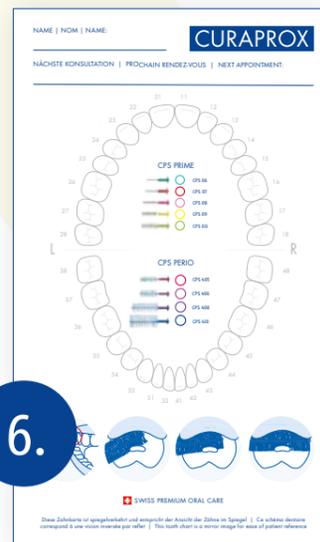


Insert the IAP into the gaps until it stops moving.


2.



3.

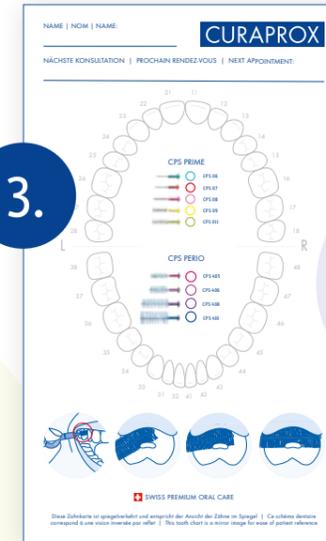

Make a note of the colour that's visible from the outside.

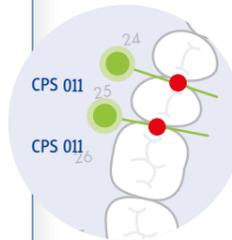

Now, take the matching interdental brush from the chairside box.


4.



The chairside box is available for CPS prime and CPS perio.

Demonstrate proper use for each interdental space. In, out, done.


5.

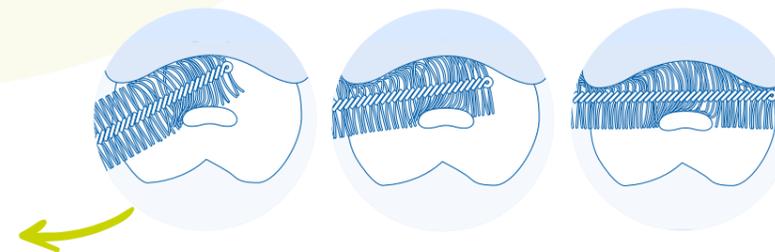


EASY TO HANDLE

All Curaprox interdental brushes can be attached to all Curaprox holders. With most holders it's as simple as: Click the old brush away. Click the new brush on. Away you go. Super convenient. Less waste. And IAP can also be clicked in without any extra hassle.

The click system

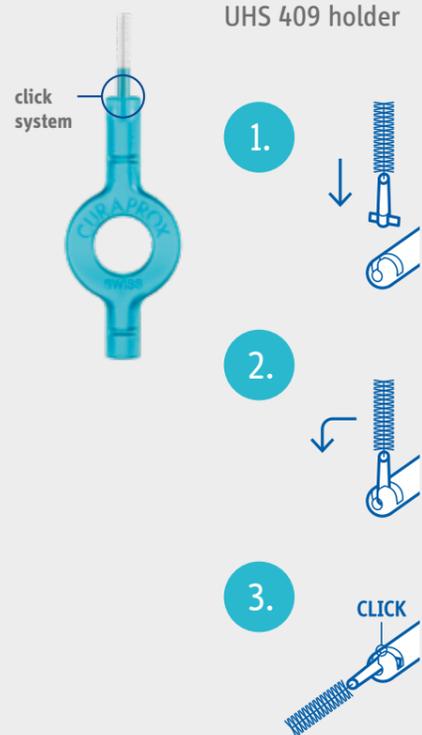

UHS 409 holder

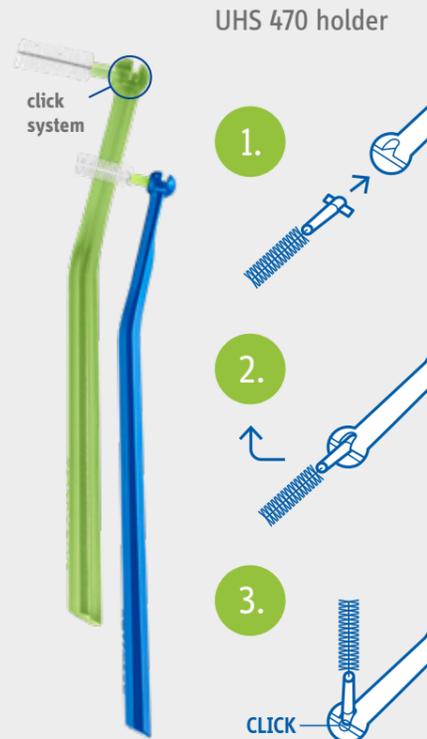

UHS 470 holder

Small, medium or large?

Made of plastic featuring the click system. Or aluminium with a twist-lock connection. Whatever your preference, all brushes fit all holders.

Less waste

Curaprox interdental brushes have a long lifespan – thanks to Cural® wire and Curen® filaments. But after 10-20 days of use, they no longer clean perfectly. Do you have to throw it all away? No. With the click system, you can simply replace the brush and reuse the holder. No stress. No hassle. Less waste.

Even more ecological

Our new refill packaging can be recycled and consists entirely of PET. It also provides even better protection for the interdental brushes, which can be clicked directly out of the packaging and clicked back in again. It's incredibly simple and hygienic.

Environmentally friendly refill packs


Our simple-to-use, hygienic new refill packaging consists entirely of PET and can be widely recycled. It also provides even better protection for the interdental brushes, which can be clicked directly out of the packaging – and just as easily clicked back in again.

Easy to measure, for optimal results


Why your patients will see results

Motivation brings success – in the short-term and long-term. The 'CPS prime', 'CPS perio', 'CPS ortho' and 'CPS implant' interdental brushes are easy and pleasant to use. In closed approximal spaces, the papillae are uninjured. That's why the 'CPS prime' can be carefully pushed beyond the papillae. With open approximal spaces, the papillae have receded. That's why the 'CPS perio' and 'CPS implant' can be pushed directly through the approximal space. Because it's so simple to use CPS interdental brushes – and because they work so well – your patients will see results in both the short-term and the long-term.


In and out over the papillae in five steps

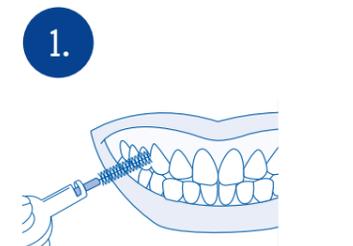


Position at an angle: Position the 'CPS interdental brush' at the approximal space at an angle to the tooth's axis.

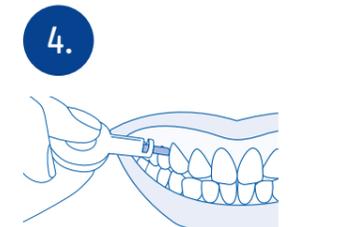


Push through: Carefully push the brush through the approximal space.

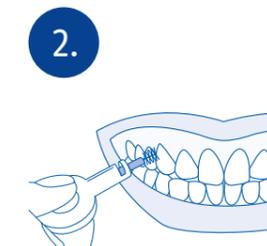

Insert to the middle: Now carefully insert under the contact point around the middle of the tooth.

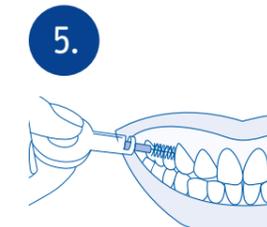

Pull out: Hold the brush at a slight angle and pull out via the papillae.

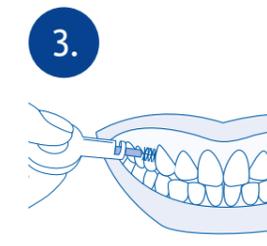


Hold perpendicular: Without moving the brush any further forward, hold the holder perpendicular to the tooth's axis.

Please note

Due to the quality of CPS interdental brushes, one movement is enough: in, out, done. Once a day is all you need.

Do not force between the teeth. If used improperly, the gums may be damaged.

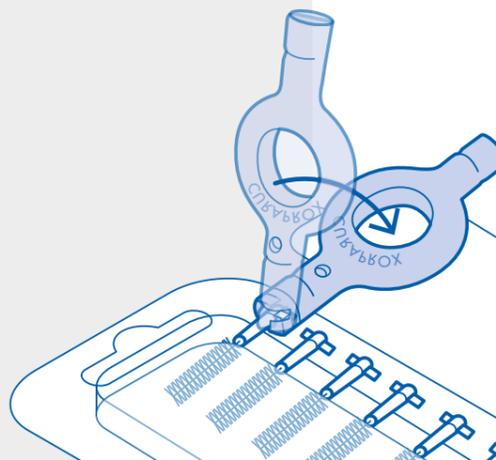

Go to curaprox.com for full instructions & video guide Professionals:

Go to itop-dental.com to learn more with the iTOP Introductory Seminar.

PRIMARY PROPHYLAXIS

The spaces in between the teeth should not be an afterthought for your patients. 'CPS prime' interdental brushes have been developed to care for healthy or only slightly injured papillae – and they've been proven to be more effective than dental floss. The ultra-fine bristles become taut to fill up the space and clean the sulcus as well as the critical areas up to the contact points effectively removing plaque and pathogens from the approximal space. This helps to prevent cavities, gingivitis and periodontitis – and even heart and systemic diseases. Patients find them pleasant to use and the risk of injury is very low. Dentists have everything they need with probes and five sizes in the chairside box.

06

CPS 06
0.6 mm* / 2.2 mm**



07

CPS 07
0.7 mm* / 2.5 mm**



08

CPS 08
0.8 mm* / 3.2 mm**



09

CPS 09
0.9 mm* / 4.0 mm**



011

CPS 011
1.1 mm* / 5.0 mm**



CPS prime

- Primary prophylaxis
- System with five sizes, IAP, IAC, colour coding
- Super-fine, extra-long, ultra-resilient filaments
- Cural® wire – extra strong, extra thin
- Colour coding
- Our click system for all holders is kinder on the planet
- Developed and produced in Switzerland


Want to find out more about interdental brushes for your perio patients? See page 62.

Interdental access probe (IAP) prime




Should I recommend dental floss or interdental brushes to my patients?

Dental floss cannot reach crucial areas.
Dental floss glides over gaps:

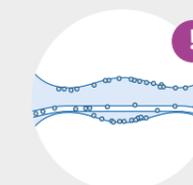

But not all interdental brushes are right for the job. You need an interdental brush in the right size and with long, ultra-fine, resilient bristles. Introducing the 'CPS prime'

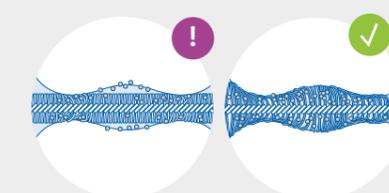

The 'CPS prime' becomes taut in the interdental space – and fills the whole approximal space, reaching into gaps and the sulcus to tackle plaque and pathogens in all critical niches. In one recent study*, participants who used interdental brushes developed significantly less plaque than participants using dental floss.

* Luz M, Klingbeil MFG, Henriques PSG, Lewgoy HR (2016): Comparison between Interdental Brush and Dental Floss for Controlling Interproximal Biofilm in Teeth and Implants. Dent Health Curr Res 2:3. doi:10.4172/2470-0886.1000119

Another study**, back in 1998, concluded interdental brushes remove plaque more effectively than dental floss.

** Christou V, Timmerman MF, Van der Velden U, Van der Weijden FA (1998): Comparison of different approaches of interdental oral hygiene: interdental brushes versus dental floss. J Periodontol. 1998 Jul;69(7):759-64. PMID: 9706852. DOI: 10.1902/jop.1998.69.7.759

* Accessibility
** Effectiveness

POWERFUL PREVENTION VERY GENTLE

In perio patients, the interdental space is open. The papillae can be seriously injured and there may be black holes. 'CPS perio' are perfect in such cases thanks to their high-performance cleaning. They're solid and resilient enough to efficiently and gently clean as secondary and tertiary prophylaxis. 'CPS perio' are also ideal for fillings, crowns, bridges and stubborn plaque. The usual 'CPS system' applies simply establish the right size using the IAP probes in your chairside box.

405

CPS 405
1.3 mm* / 5.0 mm**



406

CPS 406
1.7 mm* / 6.5 mm**



408

CPS 408
2.2 mm* / 8.0 mm**



410

CPS 410
2.8 mm* / 10.0 mm**



Interdental access probe
(IAP) perio



CPS perio


- Powerful prevention. Very gentle.
- For secondary and tertiary prophylaxis
- System with four sizes, IAP, IAC, colour coding
- Super high-performance cleaning
- Painless
- Extra resilient yet gentle filaments
- Strong wire core
- Colour coding
- Our click system for all holders is kinder on the planet
- Developed in Switzerland.
Made in Germany.

Need extra help for perio patients?
'Perio plus' gives antibacterial support.
See page 10



There are two crucial criteria for perio care: **Comfortable. And effective.**


Caring for the approximal space is particularly important for perio patients so it's important to give them the right tools

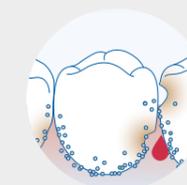

Space-filling interdental brushes are needed for secondary and tertiary prophylaxis. A certain level of rigidity is required. It's not just about being slightly harder than primary prophylaxis interdental brushes - more tension is needed.

Finding the right balance between rigid and pleasantly gentle is the key. If the filaments are too rigid, they may cause further injury. If they're not rigid enough, they won't be resilient enough to clean thoroughly.

After extensive development, Curaprox has found the perfect blend. Which means that, for the first time, perio patients and others can now benefit from the ultimate interdental brushes for open interdental spaces. Pleasant to use and extremely effective.

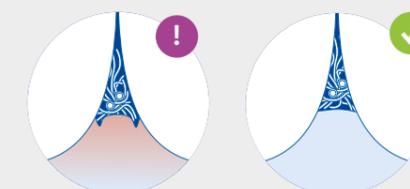

Start your patients on interdental cleaning with 'CPS perio' as soon as they've healed. You can help them to improve their quality of life. Show your patients just how easy and pleasant oral care can be with the 'CPS perio' to give them the foundation for a pleasant and successful daily practice.

* Accessibility
** Effectiveness

CARING FOR LARGER INTERDENTAL SPACES

After periodontal or implant treatment, some interdental spaces are much larger than usual. We developed 'CPS implant' – with long, fine, rounded filaments – to gently and effectively clean these spaces. It's perfect for the job. The wire core is resilient and durable and it also comes with a plastic coating to protect the surface of implants. Gentle on implants, crowns, bridges and after surgery.

505

CPS 505
>2.0mm* / 5.5mm**



508

CPS 508
>2.0mm* / 4.5-8.5mm**



CPS implant

- Perfect for cleaning large gaps
- System with five sizes
- Extra-long, extra-fine, extra-rounded filaments
- Strong wire core
- Plastic coating protects the surface of implants
- Our click system for all holders is kinder on the planet
- Developed in Switzerland. Made in Germany.

Implants make it even easier for bacteria to settle in the sulcus. The CS 708 toothbrush is the perfect solution. See more on page 45.

* Accessibility
** Effectiveness

Reach for the stars

The CPS chairside box is now available for the 'CPS implant': It contains all five sizes – ready mounted on holders.

How to use interdental brushes – the short version

Position the 'CPS interdental brush' at a slight angle in front of the interdental gap – the space between the teeth. Push it in: Just a little at first, then straighten it and push it right through. Pull it out: Very gently. Done. Move on to the next one. Go to curaprox.com for full instructions & video guide

Tip for prevention and patient retention

Get hold of our publication for patients, 'why implants need even better care than natural teeth', to help them understand why perfect cleaning and regular checks are necessary. Get one from your Curaprox sales representative. Go to www.curaprox.com

Show your patients the difference and they will understand why perfect cleaning and regular checks are necessary.

Plastic coating protects the surface of implants

Caring for implants safeguards general health

Your patients need to truly understand that artificial dentures require special care.


Part of your job as a dentist is to let your patients know that the surfaces of the implant and hard-to-reach spots must be carefully cleaned – especially the edge of the gum around the implant.

We're convinced peri-implantitis and implant losses can be avoided in many cases. When patients better understand the exact differences between implants and natural teeth, they clean more enthusiastically and more effectively.

So be sure to explain to them in detail that:

- An implant lacks the tooth fissure or periodontal fissure, which have fine blood circulation and function as a defence system
- Without this blood circulation, fewer antibodies are available to fight bacteria
- Some of these bacteria in the oral biofilm are pathogens that cause heart and systemic diseases

Regular implant care creates the conditions to prevent inflammation and peri-implantitis. Give your patients the right medical information in an easily understandable way, encourage diligent oral hygiene and regular dental check-ups with you and you can help them to enjoy implants that last for decades.

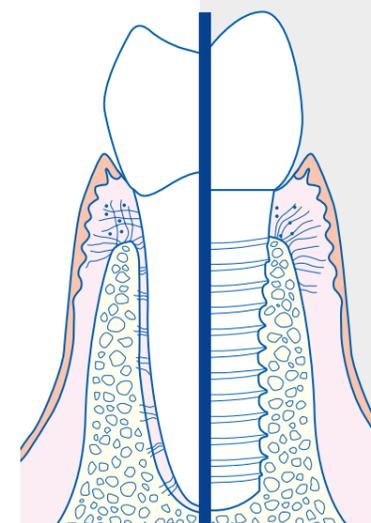

CARE FOR BRACES

Braces may mean beautiful teeth, but they also bring a whole lot of responsibilities. It's not just teeth that have to be cleaned perfectly every day – the brackets and wires do too. This is where the interdental brushes CPS 14 ortho and CPS 18 ortho can help.

14

CPS 14
1.4 mm* / 5.0 mm***



18

CPS 18
2.0 mm* / 8.0 mm***



CPS ortho


- Specially designed for ortho patients to keep spaces around their braces clean.
- Cone-shaped CPS 14 cleans orthodontic bands from the inside and outside, as well as areas between the wires and teeth and around the brackets.
- CPS 18 is ideal for cleaning the wires, simply by inserting it underneath

Implants make it even easier for bacteria to settle in the sulcus. The CS 708 toothbrush is the perfect solution. See more on page 45.



For orthodontic patients. Specially designed interdental brushes to clean wires from the inside and outside, and around brackets. Recommend the following procedure:

1. The CPS 18 ortho is ideal for cleaning the outer wires.
2. To clean the areas created by the braces, run the CPS 18 ortho under the wires.

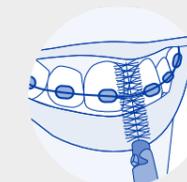

3. The best thing for cleaning inner wires is the conical CPS 14 ortho, which your patients may even find suitable for outer cleaning too.



4. It may perhaps require more than one cleaning movement.

* Accessibility
** Effectiveness

FOR PERIO PATIENTS

Periodontitis is an inflammation of the tissues that hold in place and support the teeth. If left untreated, patients risk losing a tooth. Here are some special, essential helpers to give your perio patients the tender, loving care they deserve:



Perio plus

See page 10+
'Perio plus' contains CHX together with Citrox® natural bioflavonoids – extracted from bitter oranges – to reinforce the antibacterial effect. It's effective against the causes of periodontitis and other diseases.



CPS perio

See page 62
The 'CPS perio' interdental brush is recommended for larger interdental spaces. Extremely effective in combating stubborn plaque... without causing pain. Available in four sizes. Once a day. Once in and out. Done.



CS 5460

See page 42
Brushing your teeth is pure joy. This toothbrush is so gentle that it polarises. The sulcus loves it. The biofilm hates it.



Enzycal 1450

See page 72
Gentle, mild, effective: The enzymes in Enzycal provide reinforcements... supercharging the natural power of our saliva.



Hydrosonic pro

See page 18
Professional teeth cleaning with sonic technology and innovative Swiss brush heads – for exceptional oral care. There's no hiding place for the bacteria and plaque.



CS 1006

See page 44
Tooth brushing doesn't get more precise than this. The CS 1006 – with its single tuft brush head – is the ultimate in tooth-by-tooth and gumline cleaning.

Be you

See page 74
Toothpaste available in six tastes. Six zesty flavours. Six vibrant colours. Long-lasting freshness. Gentle, natural whitening. A new dimension in oral health.





MADE WITH LOVE TO CLEAN & PROTECT


Enzycal. Gentle, mild, effective. Forget the usual harsh ingredients. Toothpaste that tastes fresh and fights tooth decay – with enzymes to boost the protective power of your saliva.

Secondly... **'Be you'**. Zesty, fresh flavours. One toothpaste. Six astonishing flavours. Choose depending on your mood. Made with joy and love, and natural ingredients.

And then there's **'Black is white'**. It's as black as night. And it's a whitening miracle. Activated carbon gently removes discoloured particles and leaves the mouth feeling fresh and clean.

SO GENTLE. SO EFFECTIVE.

You'll notice it straight away. Enzycal toothpastes hardly foam at all... because they're free from SLS – the aggressive surfactant that can unfortunately still be found in most toothpastes. It attacks the cells of the oral mucosa... and leads to painful irritation, injury and mouth ulcers. SLS, Triclosan and other nasties are not what we want for you or your patients. Enzycal cleans effectively, is gentle on the oral mucosa and can even significantly reduce mouth ulcers. We only use ingredients that work naturally with your body.... such as three enzymes that promote oral health, stimulate saliva production and reinforce the saliva's natural protective and restorative functions.

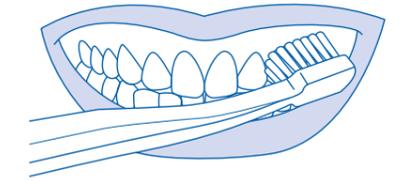

How to use: Moisten toothbrush, apply a pea-sized amount of Enzycal and you're ready to go!

- **Top tip for patients**
Recommend Enzycal to patients with dry mouth or mouth ulcers.
- **Tip: how much toothpaste should I use?**
Even for adults: a pea-sized amount of the Curaprox toothpastes Enzycal, 'Black is white' and 'Be you' is enough.



Enzycal 1450

- Suitable for age 12+
- Fresh taste
- Contains a system of enzymes
- 1,450 ppm sodium fluoride
- Low RDA



Enzycal zero

- Mint-free – compatible with homeopathic treatment
- No fluoride
- Contains a system of enzymes
- Low RDA

Enzycal

- So gentle. So effective.
- Contains a system of enzymes
- Free from SLS, triclosan, microplastics
- Low abrasion, enamel-friendly
- Perfect with manual and sonic toothbrushes
- Made in Switzerland

Enzycal is Swiss made – like our 'Be you' and 'Black is white' toothpastes. See pages 74 and 77.



Enzycal 950

- Suitable for age 2+
- Extra mild flavour
- Contains a system of enzymes
- 950 ppm sodium fluoride
- Low RDA

Enzymes – saliva's secret weapons

The saliva's protective enzymatic function can be slightly impaired – not just by alcohol and nicotine, but also by everyday life. For example, by drinking sugary, acidic drinks.

Enzycal toothpastes support the saliva's natural protective system. Enzycal contains – at a physiological level – just those enzymes that promote oral health. This enables saliva to perform its vital functions in the mouth to the fullest.

These enzymes are amyloglucosidase, glucose oxidase and lactoperoxidase. In the body they fight bacteria, fungi and viruses, prevent dental build-up and break down harmful sugar before micro-organisms can turn it into acid.

Here's how it works: Amyloglucosidase turns polyglucan, the matrix of plaque, into glucose. Now glucose oxidase comes in, turning glucose into hydrogen peroxide, which competes with plaque bacteria – reducing their ability to reproduce. Lactoperoxidase then ensures that hydrogen peroxide from the glucose transformation reacts with thiocyanate to create hypothiocyanite – an antibacterial substance naturally found in saliva.

The result: The growth of bacteria is prevented and the natural balance of oral microflora is restored. Less plaque, less tartar, better oral health.

GENTLE EVERYDAY WHITENING

Gentle care for your teeth and gums

'Be you' protects the enamel with a low RDA. Hydroxyapatite closes openings that expose dentine – and relieves the pain of sensitive teeth. Glucose oxidase encourages the flow of saliva with natural enzymes and strengthens your mouth's natural defences. No SLS means no dry mouth – and fewer ulcers.

Happy and healthy every single day

Daily support for a healthy oral microflora and daily protection against caries. Fluoride, xylitol, provitamin B5 and herbal extracts fight bacteria and help to stop inflammation. And 'Be you' comes in six zesty, fresh flavours. Pure joy every day!

Naturally whiter and brighter, day by day

This is gentle, natural, whitening – enzymatically... so your teeth become a little whiter and brighter every time they brush. Hydroxyapatite smoothes the teeth, making it difficult for new discolouration particles to take hold. No harmful bleaching. No grinding down. And a filter adsorbs yellow light to make your teeth appear even whiter. Just 'Be you'... and smile for the camera!

Be you

- Toothpaste in six tastes
- Gentle, natural whitening
- With glucose oxidase, hydroxyapatite, xylitol
- Beautiful herbs including echinacea, bitter orange, devil's claw and Indian pennywort
- Six zesty flavours
- 950 ppm sodium fluoride
- Low RDA
- Vegan
- Free from SLS, Triclosan and microplastics
- Made in Switzerland

Too colourful for you? See page 77 for 'Black is white' toothpaste with activated charcoal.



Curaprox Swiss premium oral care.

The 'Be you' formula:

- It starts with fluoride. 950ppm.
- + Hydroxyapatite, meanwhile, smoothes and brightens the enamel, closes openings that expose dentine and thus reduces the sensitivity of exposed tooth necks.
- + Birch sugar – also known as xylitol – gives 'Be you' a touch of tasty sweetness, strengthens your teeth and fights bacteria.
- + And glucose oxidase, a natural enzyme found in saliva, gently whitens the teeth...
- + ... and also takes care of your tooth enamel and gums because it activates the natural lactoperoxidase system of the saliva.
- + All of which gives a gentle yet powerful antibacterial, antiviral and antifungal effect.
- + The beautiful herbs and natural ingredients of 'Be you' – including extracts of echinacea, bitter orange, devil's claw and Indian pennywort – nourish the gums and the oral mucosa and help to stop inflammation.
- + And provitamin B5 keeps the mucosa and gums moist and helps them to build new cells.
- + 'Be you' is enamel friendly thanks to a low RDA
- + You won't find SLS, Triclosan, bleaches or plastic microbeads...
- + ... but you will find pure joy!

Gentle, natural, whitening – enzymatically

The enzyme glucose oxidase reacts with the sugar in the mouth to create a low concentration of hydrogen peroxide in the saliva that gently whitens by removing stains and discolouration. And because Hydroxyapatite smoothes the enamel, it makes it difficult for new discolouration particles to adhere to the teeth. Your tooth enamel also becomes whiter as it adsorbs Hydroxyapatite. Now you can enjoy gentle, natural, long-lasting whitening.

The 'Be you' difference

Stay fresh for hours

Feel fresh and feel confident... long after brushing. Power bubbles in 'Be you' give you a refreshing and long-lasting blast of menthol. No plastic microbeads!

Some crazy science

- A filter in 'Be you' adsorbs yellow light on the tooth's surface
- Yellow light is not reflected
- The tooth appears even whiter!

Strictly no nasties

Harsh ingredients are our enemy. For example, SLS can cause irritation and mouth ulcers. At Curaprox, we use gentle, harmless-but-effective ingredients only. No SLS, no Triclosan, no bleaches, no plastic microbeads.

How your patients can improve the whitening effect even more!

- + Don't rinse after brushing
- + Wait for an hour until eating or drinking
- + Definitely don't smoke
- + Brush just before bedtime

How to use: Moisten toothbrush, apply a pea-sized amount of 'Be you'. Brush. Feel fresh!

Top tip

Recommend 'Be you' to your patients to show them that oral health is anything but boring.

BE YOU



GIN TONIC + PERSIMMON



Brushing the sophisticated way: Sweet persimmon fruit. A twist of juniper berry. All the fun of a gin & tonic without the alcohol. Ready to embrace the challenges of the day.



PEACH + APRICOT



'Wow' for your teeth: What makes you happy? What makes you say 'wow'? Subtle notes of peach and apricot. A touch of refreshing mint. Wow. Pure happiness.



GRAPEFRUIT + BERGAMOT



Toothpaste to help you shine: Juicy grapefruit, invigorating bergamot – wake up your mouth, wake up your brain. It's time to shine like a star.



WATERMELON



Treat your teeth: Sweet, icy fresh watermelon. Refreshing menthol. Feel the anticipation... like a little kid looking forward to a trip to the candy store. Love life.



BLACKBERRY + LICORICE



Sumptuous luxury in a toothpaste: Dreamy blackberries meet luxurious liquorice – first thing in the morning... and last thing at night. Dream wonderful dreams.



APPLE + ALOE



Make brushing extraordinary: We explored flavour combinations. And we discovered apple, aloe vera and fresh mint. Extraordinary. Explore the unknown... you never know what you might find.

BLACK IS WHITE

ACTIVATED CARBON A SECRET SUPERPOWER FOR BRIGHTER, WHITER TEETH

It's a charcoal toothpaste. It's as black as night. But it gently – yet powerfully – whitens... without harsh, abrasive ingredients. Beautiful, whiter teeth give your patients new confidence. But normal whitening toothpastes are harmful. Curaprox – the oral health experts – give you the best of both worlds with 'Black is white' toothpaste. This is supercharged whitening – the gentle, healthy way.



Black is white

- Gentle, effective whitening...
- ... without harsh painful chemicals
- Activated carbon naturally adsorbs discolouration particles
- Enzymes support the saliva's natural functions
- Hydroxyapatite smoothes and brightens the enamel and reduces the sensitivity of exposed tooth necks
- Physical blue filter makes teeth look whiter
- 950 ppm sodium fluoride
- Free from SLS, Triclosan, microplastics
- Suitable for sonic toothbrushes
- Made in Switzerland

Want more colour in your brushing?
More atmosphere? See page 74.

Black is white
Fresh lime-mint

Activated carbon...
... is completely natural. Because it's porous, it simply adsorbs discoloured particles to give your patients high-performance whitening. And it leaves the mouth feeling fresh and clean. Gentle yet powerful, natural, effective whitening and oral care to give your patients something to smile about.

WHITENING

WE LIKE WHITE TEETH – BUT WE LOVE HEALTHY TEETH

Whitening should never damage your gums or oral mucosa. Oral care should never be harsh or abrasive. Your patients need to protect their teeth enamel. If it starts to erode, it can lead to sensitivity and pain. This is gentle, natural, healthy whitening with outstanding results:

‘Black is white’ Hydrosonic

See page 22
Dentistry meets design. Sonic technology and effective tooth cleaning with style. With up to 84,000 brush movements per minute, the Hydrosonic ‘Black is white’ can easily handle biofilm and discolourations – and it’s super gentle. The brush heads are also black. The activated carbon particles adsorb stains without damaging your enamel.

‘Black is white’ manual Toothbrush

8,760 ultra-fine Curen® carbon filaments with activated carbon particles gently remove discolourations without damaging enamel. Gentle on teeth and gums. Relentless on plaque.



Travel set

See page 47
Whitening on the go. Make sure you always have your favourite ‘Be you’ toothpaste by your side with the Curaprox Travel set. It also comes with a practical Curaprox CS 5460 travel brush and a prime interdental brush. Go!



‘Be you’

See page 74
The natural enzyme glucose oxidase converts sugar into hydrogen peroxide and then oxygen. This is natural, gentle, sustainable whitening.



‘Black is white’ chewing gum

Is chewing gum good for you? ‘Black is white’ enzymatic whitening chewing gum is. But not just because it gives you lemon-minty-fresh breath. It actually starts a flow of saliva. Chew your way to a healthier mouth and whiter teeth.

Top tip: Your patients may be shocked to learn that other whitening gums contain ingredients such as marble powder or similar abrasive materials. You can recommend a whitening chewing gum that’s actually good for them.



CURADEN ACADEMY – THE PLACE WHERE LEARNING COMES TO LIFE



Get involved:
CuradenAcademy.com

The future of dentistry depends on innovative treatment options, but also on how we pass on knowledge and skills to patients to keep their own smiles healthy for life. This is the reason Curaden Academy exists – to help dental professionals acquire the knowledge and skills necessary to be the best dentists and patient educators of the future. Our courses started more than fifteen years ago with our most famous brand of courses, iTOP (individually trained oral prophylaxis) which was developed by Jiří Sedelmayer, a dentist and clinical teacher that worked in Hamburg University for over two decades. iTOP courses and their Hands-On concept called Touch2Teach have redefined the approach to prevention, combining the latest research on oral diseases, biofilm control techniques, and behaviour change knowledge to create a leading course in oral health prevention. That, as well as the significantly increased demand for individual training and coaching, has resulted in more than 30,000 trained dental professionals in over 70 countries worldwide.

We take a different approach to learning. Whether we organise a surf camp for dental students where we combine studying and pleasure, or we create an interactive webinar with world-class leaders in their respective areas, Curaden Academy courses are there to provide the best learning experience that will benefit you, your team, your practice and your patient.

Today, Curaden Academy offers a wide range of courses, both live and online, that accommodate the needs of dental professionals worldwide. We collaborate with more than 30 international dental universities, several dental associations and we are proud to work with dental experts in every corner of the world.

Join us on <https://CuradenAcademy.com> and check out our latest courses and events!



'We need to understand each patient's uniqueness in order to select an individual education strategy. That's the key to effective oral health prevention.'

Mário Rui Araújo PhD, M. Psych, BSDH,
iTOP Lecturer from Portugal

'With iTOP, I acquired exclusive knowledge and skills that I pass on directly to my patients – along with motivation and enthusiasm. My patients know it is worthwhile and also worth the money, they love it.'



Dr Rolf Kufus, works with the
iTOP philosophy, Switzerland



'iTOP – The 'lightbulb moment' in preventive dentistry.'

Elna Van der Ham, RDH,
iTOP Lecturer from South Africa



'Good dentistry focuses on prevention – not just cure. Having natural teeth beats implant, crown or bridge.'

Dr. Carlos Molina, Periodontist,
iTOP Lecturer from Ecuador

'Curaden Academy courses made me see that with good prevention, there would be no need for restorative treatment.'



Prof. Dr. Hugo Lewgoy,
Dental Surgeon,
iTOP Lecturer from Brazil

Even more training:
CuradenAcademy.com

CURADEN CORP.: 'BETTER HEALTH FOR YOU'

In 1954, Hans Breitschmid asked the question:
'How can I help dentists to do their jobs better?'

The question was asked so many times over the years and decades that it eventually became Curaden Corp. mission – and the reason for its success. More and more answers were found. Answers that anticipated the future. For this we can thank the talent of Hans Breitschmid – and later his son Ueli – for listening to their customers and dental professionals... and being the first to anticipate changes and deliver the best new solutions and processes.

For example: The world's first treatment unit for ergonomic work on recumbent patients. Probably the world's first interdental brushes for primary prophylaxis. The first CHX mouthwash containing Citrox®. Perhaps the first toothpaste to contain gold. And the first courses to teach 30,000 dental professionals worldwide unique methodology in oral healthcare.

Hans Breitschmid's question has expanded to become: 'What is a dentist's job?'

The whole dental industry struggles to answer this question. Curaden Corp. answer? A dentist is a preventative practitioner who helps the patient to achieve better health. Health that begins, but does not end, in the mouth – the domain of dental medicine. For example, it's crucial to balance the microbiome in the mouth to prevent periodontitis bacteria from entering the bloodstream... thus preventing heart attacks and strokes. The mouth is the gateway to health. The dentist is the gatekeeper.

Then we have society's inherent need for beauty. In some cases, aesthetic dentistry is appropriate. But the best oral care makes the biggest difference for both health and beauty... especially when monitored by an exceptional dentist who always finds the time to instruct, coach, motivate and recommend the best products. That means keeping the mouth and the teeth healthy for life – ideally completely free of cavities, implants, bridges and prostheses. A healthy mouth for life.

To help our dentists and customers achieve this goal, we offer exceptional products from Curaprox and exceptional training with iTOP and our BOB App.

Curaden Corp. is a Swiss family company based in Kriens, Luzern. We work closely with dental professionals to promote our Curaprox products... all invented and developed in Switzerland. Attractive, effective products that are a pleasure to use. They are available in practice from dentists and dental hygienists, pharmacists and chemists as well as online at www.curaprox.com.

CURAPROX

curaden
better health for you

Curaden AG
Amlehnstrasse 22
CH-6010 Kriens
T +41 (0)41 319 45 50
F +41 (0)41 319 45 90
info@curaden.ch

www.curaprox.com



 SWISS PREMIUM ORAL CARE