

# CURAPROX TOOTHBRUSHES



**CS 12460 velvet**  
 12,460 special Curen® filaments guarantee the smoothest Curaprox experience of all. For strong beautiful teeth and healthy gums... without enamel damage, without receding gums.



**CS 5460 ultra soft**  
 The top-selling Curaprox toothbrush. Gentle on the gums and extremely efficient at cleaning. Trusted and loved around the world.



**CS 3960 super soft**  
 With 3,960 Curen® filaments: not ultra soft, yet still super soft.



**CS 1560 soft**  
 For people used to hard toothbrushes. The ideal transition to the softness of Curaprox.



**CS Smart**  
 For children aged 4-12 and adults who prefer a smaller brush head. With the finest 7,600 Curen® filaments.



**CS ortho ultra soft**  
 Special ortho brush: special brace-friendly design.



**CS 1009 single**  
 Special ortho brush: extra long bristles for perfect plaque control even in hard-to-reach areas.



**CS 1006 single**  
 Tooth by tooth. Brushing doesn't get any better than this. The perfect toothbrush for expert users.




**CS surgical mega soft**  
 The mega soft brush for extremely sensitive gums especially after operations, inflammation, and radiotherapy in the head area. Ideal in combination with chlorhexidine.



**CS 708 implant | ortho**  
 Super compact brush head for perfect plaque removal from implants, fixed prosthetics and braces.




**CS kids ultra soft**  
 For children aged 4-12. With the finest 5,500 Curen® filaments.

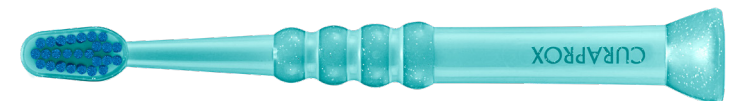


**Baby toothbrush**  
 For children aged 0-4. Gentle on the gums: thanks to a rubberised brush head and soft 4,260 Curen® filaments.

# CURAPROX

## THE SECRET BEHIND CURAPROX TOOTHBRUSHES

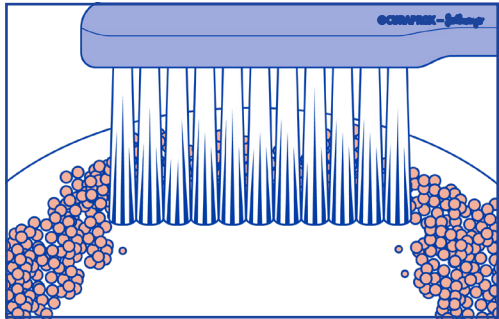

**Curaprox toothbrushes** remove and disrupt plaque thoroughly - thanks to their exceptionally dense bristles.

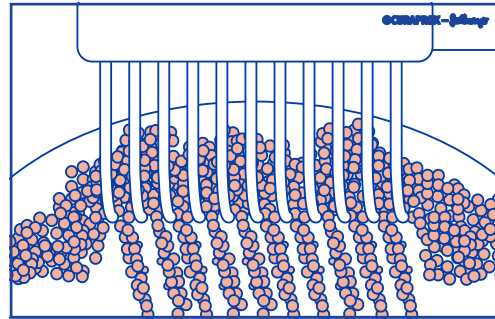


**Normal brush with few bristles:** Plaque remains on the teeth between the individual bristles.



- Finest Curen<sup>®</sup> filaments  
→ gentle on your gums
- Tightly packed Curen<sup>®</sup> filaments  
→ effective against plaque

## The right technique

### Brushing angle

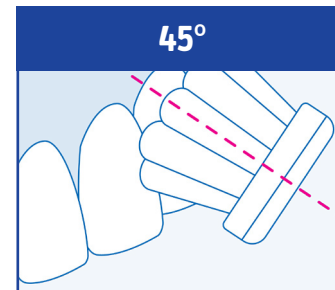

### Pressure

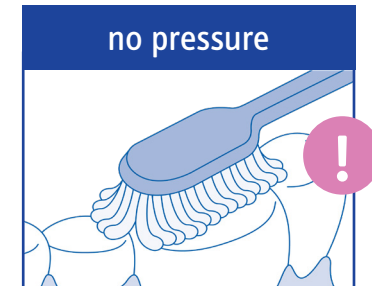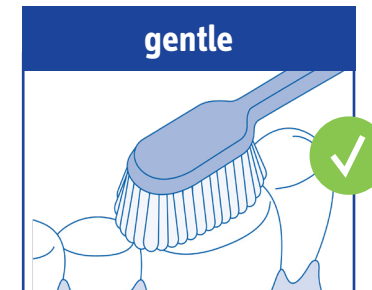

Maintain a 45° angle: half on the gum, half on the tooth, and brush in small circles without applying any pressure.